

Phase Transitions:
Metal–Insulator Transition in VO_2 and
Ferroelectric Transition in BaTiO_3

Advanced Laboratory
Department of Physics & Astronomy
Stony Brook University

Write-up version of September 2026

Sections on ferroelectricity in BaTiO_3 by Laszlo Mihaly and, in part, Michael Gurvitch; sections on phase-transition theory and the metal–insulator transition in VO_2 by Michael Gurvitch.

Contents

Pre-experiment and mid-experiment evaluation sheet	4
1 How to prepare for the measurement	5
2 Goals of the experiment	5
3 Experimental set-up	6
4 Principle of the measurement and an overview of data analysis	6
5 Looking ahead: what should and should not be included in your Final Report	8
6 Measurement techniques	9
6.1 How to use lock-in to measure resistance?	10
6.2 How to use the lock-in amplifier to measure the capacitance?	13
6.3 Lock-in amplifier settings	14
7 The rest of the experimental apparatus in detail	17
7.1 Sample chamber	17
7.2 Samples and sample preparation	17
7.3 Temperature control	20
7.4 Temperature measurement	22
8 Suggested protocol for a measurement	23
9 Advice for evaluating the capacitance data	23
A Elements of AC circuit theory relevant to our measurements using the lock-in amplifier	25
A.1 Measuring capacitance	26
A.2 Alternative derivation using complex numbers	28
B The four-probe method: advantages and limitations	30
B.1 How much voltage-contact resistance can be tolerated?	32
C Data taking with PyFERRO	35
D The EG&G 5302 lock-in	39
E Physics of phase transitions	40
E.1 Introduction to phase transitions	40
E.2 Classification of phase transitions	41
E.3 Thermodynamic approach to phase transitions	43
E.3.1 First-order transitions from the thermodynamic point of view; mixed-phase states	43
E.3.2 Second-order transitions from the thermodynamic point of view; elements of Landau theory	45
E.4 Description of the phase transitions and materials you will study	48
E.4.1 Metals, semiconductors and insulators, and their conductivity	48
E.4.2 Metal-insulator transition	51
E.4.3 Essential physics of one type of MIT: the Mott transition	51
E.4.4 MIT in VO_2	52

E.4.5	Development of a mixed state and hysteresis in VO_2	54
E.5	BaTiO_3 and ferroelectricity	55
E.6	Paraelectric to ferroelectric transition; Curie–Weiss law	57
E.7	Derivation of C_{CW} in terms of a microscopic dipole moment	58
References		62

Pre-experiment and mid-experiment evaluation sheet

You must meet with a staff member for the experiment you are doing before beginning an experiment and before the lab period indicated in the middle of the experiment for an evaluation of your preparation and your progress to date. This form must be signed and dated by the staff member to indicate you have done the required work to date. It must be attached as the first page of your report. You will lose one grade point from your report grade for each section of this form that is not properly signed on time. If this is your first experiment you should also be prepared to demonstrate you have read and understood the information on log books in the course notes.

Student name _____

Partner name _____

Pre-lab discussion

Be prepared to demonstrate you have read the background material for this experiment (most of this material is contained in this write-up, but also in cited references; it includes theory described in Appendix A). Discuss the following questions with the staff:

1. Describe what we mean by phase transitions; give some examples; what is the difference between phase transitions of the first and second kind? Experimentally, in this lab, how can we tell which transition it is?
2. Give a general description of a metal–insulator transition; what happens to the electronic spectrum of a solid when this transition takes place? What is the most direct manifestation of this transition, and how can it be detected? Is this transition taking place all at once in the whole sample, or not (explain)?
3. Explain what ferroelectricity is and what the mechanism for it in BaTiO_3 is.
4. What does a PID controller do? Be specific about its principle of operation.
5. What is the general principle of operation of a lock-in amplifier?

Completed before first period for the experiment? Yes / No

Staff comments _____

Staff signature _____ Date _____

1 How to prepare for the measurement

This write-up contains a lot of material. You should study parts of it before coming to the lab, and read and study the entire text throughout the course of this experiment.

Note. A big chunk of this manual is in the Appendices, so do not treat them as secondary-importance material.

You should also supplement this write-up by looking up and reading some of the books/articles/websites from the list of references. We propose that you start in the following order:

First read Sections 2–5, which will give you an overview of this lab. Then spend some time learning how to use the lock-in amplifier. Your first task will be to master preliminary practice measurements of resistance and capacitance. During the first two periods of the lab work, concentrate on **Section 6** of this write-up entitled “Measurement techniques” and study **Appendices A and B**. Read the equipment manual for the lock-in amplifier and Ref. [2]; study the principle of its operation. Simultaneously, study the instrument hands-on: connect a circuit and first measure a known resistor, and then a known capacitor, using the lock-in.

Once you get familiarized with using the lock-in, read about the rest of the experimental set-up, which is described in Sections 7–9. Then you should start reading more of a physics theory, which is to be found mainly in a large Appendix E. The physics of phase transitions is not trivial, and it will take time to “digest” and understand it.

From the beginning, keep good records in your notebook. All electrical circuits that you will build, simple calculations which accompany construction of these circuits, such as estimates of the measuring current, choice of a limiting resistor, calculations of resistance and capacitance, sketches of samples or experimental apparatus, details of sample preparation and geometry, all notes concerning physics, and all your calculations of physical quantities should be recorded in the notebook. Short comments you may write should be understandable, in case you or someone else will later examine your records. The actual experimental data will be stored in the computer and need not be in the notebook.

2 Goals of the experiment

- Understand experimental techniques used in this lab, including the principle and practical operation of a lock-in amplifier when measuring resistance and capacitance, the four-probe method, temperature control.
- Learn to handle and prepare samples; learn how to quickly check samples before they are used in lengthy measurements.
- Study the metal–insulator transition in VO_2 .
- Study the ferroelectric transition in BaTiO_3 .
- Learn basic physics associated with phase transitions and apply it to interpreting specific experiments at hand.
- We do not just do measurements for the sake of measurements. An experiment should lead to some conclusions and it should culminate and result in learning of some “new” physics (new to us in this lab; perhaps not new to the physics community). Therefore it is very important to use your data to extract some meaningful physical results: namely, estimates of microscopic or theoretical quantities which enhance our quantitative understanding of the studied phenomena. You will do so for both phase transitions studied in this lab.

3 Experimental set-up

- Sample chamber with four electric leads.
- Temperature controller.
- Lock-in amplifier (SRS SR830) to measure the samples.
- HP 34401A multimeter (not needed for the standard measurement).
- Computer (username `LabStudent`, password `PhaseTransition`).
- Decade resistance box, additional resistors, capacitors, BNC cables etc.

The Omega CN08D3-SR-2U-CM-AC temperature controller controls the heater and can be programmed from its front panel, connected to the computer via RS-485. The lock-in amplifier is connected to the computer by the IEEE 488 (GPIB) interface and PyFERRO displays and records all measurements in real time.

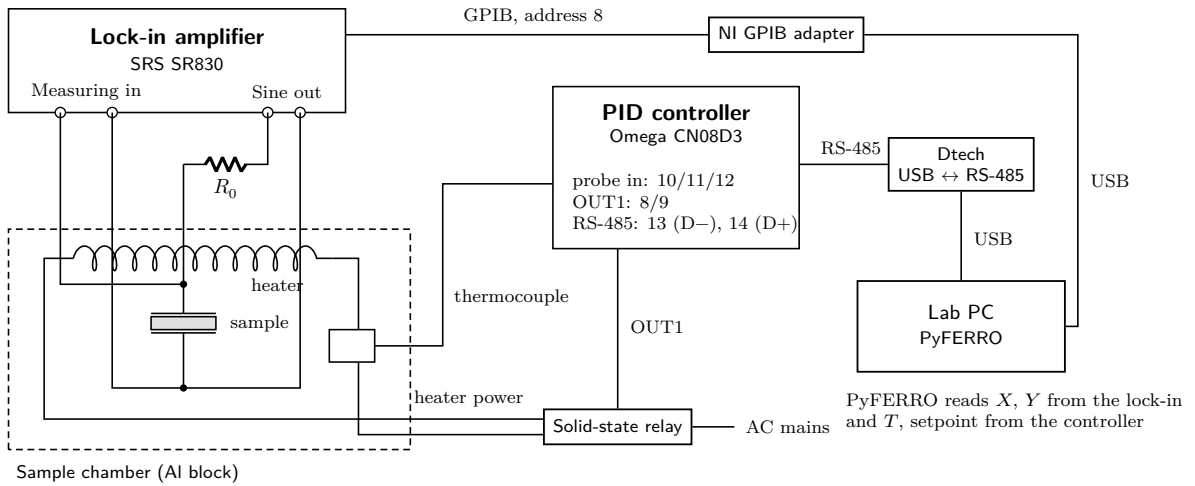


Figure 1: Instrument connections for the temperature-dependent measurement.

4 Principle of the measurement and an overview of data analysis

Details on these subjects will be found in the rest of this write-up; and you know what they say: *the devil is in the details*. Here is the **brief summary**:

You will be studying two different phase transitions, both taking place as a function of temperature: a metal–insulator transition in vanadium dioxide (VO_2) and a ferroelectric transition in barium titanate (BaTiO_3). Read **Appendix E** where the physics of phase transitions is reviewed and detailed.

The sample is placed in a temperature-controlled chamber. The temperature is slowly changed, first rising, and then cooling down, while either resistance R (in VO_2) or capacitance C (with BaTiO_3 serving as dielectric) is measured and data is recorded. The temperature will change between room temperature and about 100°C in the case of VO_2 , and between room temperature and about 150°C in the case of BaTiO_3 . The measurement of R and C is performed in an AC mode using a lock-in amplifier, which is described below in Section 6 and in Appendix A.

What will you be able to learn from these measurements? In the case of the resistive VO_2 data, you will be able to make conclusions concerning the nature and order of a phase transition,

including estimates of the metal–insulator transition temperature, of the percolation transition temperature, of the coercive temperature. You will also be able to analyze the data to extract the gap of VO_2 in a semiconducting phase and to estimate absolute resistivities of the two phases.

In the case of a capacitance measurement, BaTiO_3 will be the dielectric inside a capacitor. In a parallel-plate capacitor,

$$C = \frac{\varepsilon A}{d} = \frac{\varepsilon_0 k A}{d}, \quad (1)$$

where A is the area of the capacitor, d is the distance between the plates, ε is the dielectric constant, ε_0 is a constant called the permittivity of free space, and k is the relative dielectric constant. We will be interested in finding k as a function of temperature. It will not be possible to exactly measure A and d in the type of samples you will have, but they can be roughly estimated. So, by measuring C (which is measured quite precisely), you will be able to estimate k to within about a factor of two accuracy, the uncertainty coming almost entirely from a limited ability here to evaluate capacitor geometry. At the same time, you can certainly follow how k will change with temperature with the full precision of the capacitance measurement. This T -dependence will pinpoint the temperature of a phase transition, T_C . Combined, these data will allow a calculation of microscopic parameters of a ferroelectric transition as outlined in Appendix E. At the end, you will be able to estimate an interesting microscopic quantity: the very small distance by which atoms shift in the lattice of BaTiO_3 in a ferroelectric transition. The ability to estimate such a quantity is what makes this experiment especially interesting.

List of physical quantities which you will be expected to extract, calculate and estimate from your data

(Some of these can be measured rather precisely; others only roughly estimated.)

In VO_2 :

- The semiconducting gap E_g
- The resistance ratio between semiconducting and metallic phases
- Absolute resistivities in semiconducting and metallic phases
- The value of T_C
- The coercive temperature T^*
- The temperature of a percolation transition T_M

In BaTiO_3 :

- T_C
- Relative dielectric constant k vs. T , from about 20°C to 150°C
- Curie–Weiss law at $T > T_C$ and Curie–Weiss constant C_{CW}
- The value of a microscopic dipole moment p
- The effective atomic displacement x responsible for the appearance of p
- The last two quantities can be also calculated for BaTiO_3 single crystals using the literature value of C_{CW}

5 Looking ahead: what should and should not be included in your Final Report

Your performance will be judged primarily from the contents of your Final Report. In it there will be several different (but related) parts:

- A physics introduction concerning phase transitions and materials which you studied and a physics analysis of your results at the end. To get these right, you have to study [Appendix E](#) of this write-up and related literature.
- A description of experimental set-up (including electrical circuit diagrams) and experimental techniques that you used (lock-in and how it was used; four-probe method), including significant details of how you obtained your data. For example, if you are measuring some voltage, but your “final product” is capacitance, you will need to explain in some detail how you proceeded from that voltage to knowing that capacitance. This may include phasor diagrams and some formula derivations given in [Section 6](#) and [Appendix A](#).
- Initially you will be measuring *known* resistances and capacitances just to learn how to connect the circuits and how to use the lock-in amplifier. It is very important to perform such calibrations or baseline measurements in order to be sure that later you will be measuring your samples correctly. While a detailed description of such calibrations may be perfectly appropriate for your mid-lab report, it should not take center-stage in the Final Report, except as a short mention of the fact that these calibrations were done. You may want to report the uncertainty (i.e. % of relative deviation) in the measurement of a known R and C (of course, assuming that the values of these R and C were known very accurately in the first place). But detailed data on such baseline measurements would be out of place; the focus of the Final Report should be on measuring real samples which exhibit phase transitions.
- You will present results of sample measurements, your data, in the form of graphs, possibly some tables (but do not list all of the data points in long tables — this is neither necessary nor desirable in the Final Report). If you are measuring resistance R and capacitance C as a function of temperature T , your data may consist of graphs of $R(T)$ and $C(T)$.

Let us discuss this in a little more detail. There exists a notion of *raw data*, the data you directly obtain in your measurement, directly read off the instruments. And then there is data which was to some degree analyzed and organized. For example, you may be measuring separate values of voltage V and current I , and later you may be combining the two when obtaining resistance $R = V/I$. But your interest — and the interest of the reader of your report — is in that resistance R , and not in V or I . Thus one should use common sense when reporting data: probably in this example, separate values of V and I should not be reported (they can remain in your electronic data storage, to be examined if necessary), and R should be reported. Likewise, if you are measuring a voltage signal which is being converted to a capacitance, do not report those voltages; report capacitance directly. You may however mention the range of the voltage signals that you measured in order to tell your reader something meaningful about the quality of your data.

- The data should be presented with uncertainties, which not only should be included, but accompanied with a brief explanation of how they were obtained or estimated. If a calculation is involved, the error must be propagated from the raw data to the calculated quantity. If curve fitting is involved, the uncertainty will come out of the fitting routine. You need also to point out which uncertainties are more significant, and which can be ignored; in this lab this last comment will be particularly relevant, as some of the uncertainties will far exceed others.

- There will be data analysis: corrections which you introduced to your data, extraction of physical parameters of interest. And at the end there will be some conclusions based on your results. You will need to study this write-up and references in order to “connect the dots”, to make your experiments meaningful. When applicable, you will compare your results to the published results. Finding relevant publications by doing literature searches, reading them, and extracting from these publications important parameters, is part of a skill of an experimental physicist.

6 Measurement techniques

The main instrument you will be using in your measurements in this lab is called a *lock-in amplifier* [2]. This wonderful instrument was invented in the 1950’s by a physicist Robert Dicke; it is used extensively for a variety of measurements in physics experiments and elsewhere.



Figure 2: Robert Dicke (1916–1997), a famous American physicist who made important contributions to the fields of astrophysics, atomic physics, cosmology and gravity. He also invented the lock-in amplifier in the 1950s, with the first commercial model developed in 1962.

The advantage a lock-in gives to an experimentalist is in that it rejects a good deal of noise, increasing the signal-to-noise ratio and allowing the use of small currents which will not damage even the most delicate of samples. Additionally, by virtue of always measuring in an AC mode, it automatically averages the readout voltage, getting rid of the so-called thermal emf offsets and DC level drifts which are endemic with DC measurements, even when we are using the best and most accurate DC voltmeters.

While the lock-in amplifier was invented to measure weak signals in the presence of large noise, in our measurements we will not be pushing the limits of its capacity in that regard. Instead, we will be using its ability to measure signals of different specified phases, as explained below and in Appendix A. This lab gives you an opportunity to learn about this important instrument.

The name “lock-in” refers to locking to a frequency and to a phase of an AC signal, and measuring only that part of a signal which is in a specified phase relationship to the reference. This is a version of the *phase-sensitive detection* technique. The phase can be continually varied and also changed in 90° steps.

You should read about lock-in operation in our lock-in manual [1]. Here and in Appendix A we give a brief general introduction to its operating principle, and discuss in more detail the two types of measurements we will be doing in this lab.

Review AC circuits first. If you have not been dealing with AC circuits lately, this would be a good time to thoroughly review them: read about AC currents in general, and the phasor method in particular, in your introductory physics course textbook; some of this is also discussed in Ref. [3]. Study this section and Appendix A; it is essential that you understand this topic in order to do this lab.

The principle of lock-in measurement is as follows: the lock-in contains its own AC voltage source called *oscillator output*, in which you can adjust the amplitude, and frequency (adjustable in a wide range from a fraction of 1 Hz to 100 kHz). The frequency, once you set it, will be kept with great precision and stability by a lock-in. The readout is tuned to read only voltage at this frequency (actually, at a double frequency, but this is immaterial for our general discussion), all the other frequencies being filtered out. Moreover, it will *look* into a specified *phase* of this signal. The exact meaning of this last statement will be explained below and in Appendix A.

The *measuring input* (actually, there are two of them, A and B, and you can also hook the instrument up to measure the difference $A - B$) is an AC voltmeter, which will measure rms values, and where the sensitivity of a readout can be changed according to experimental needs, from several volts to a few millivolts full scale. It can be connected across different parts of a circuit: for example, it can be connected across R or across C . Let us now discuss specific measuring arrangements: one for measuring a resistance, another for measuring a capacitance.

6.1 How to use lock-in to measure resistance?

Suppose we want to measure a resistance, which may be changing in the course of a measurement, for example, as a function of temperature. This will be the case in our VO_2 measurement in this lab. We may connect the lock-in as shown in Figure 3.

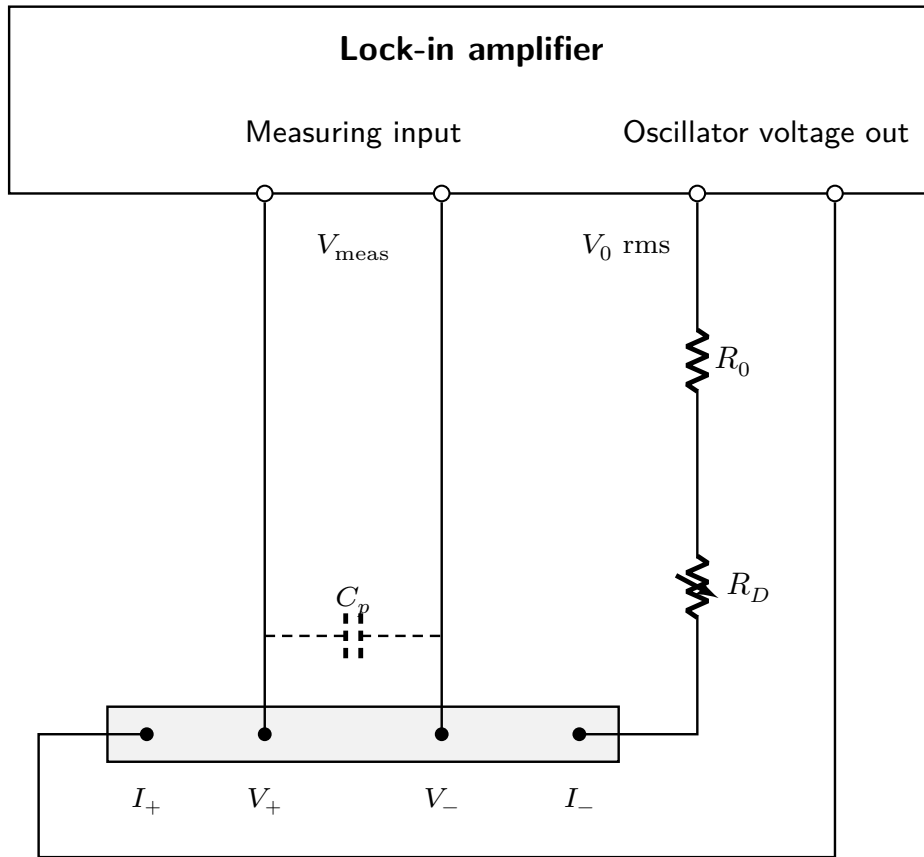


Figure 3: Connecting the lock-in amplifier for a four-probe resistance measurement. $R_0 = 50 \Omega$ is the oscillator's internal output resistance, R_D the limiting decade resistor, and C_p the parasitic capacitance of the cables. The current enters and leaves through the outer contacts $I_{\pm}$; the voltage is measured across the inner contacts $V_{\pm}$.

The label under the lock-in's SINE OUT connector shows that the oscillator circuit has internal resistance $R_0 = 50 \Omega$, which is really a part of the lock-in, but, for clarity, is shown on

the diagram above outside of a lock-in block. It is connected in series with the limiting decade resistance R_D and the sample. The measuring input of the lock-in is connected across the inner contacts of the sample.

Let us first discuss the current. It is alternating, AC current; we will be always talking about rms values of all voltages and currents in this discussion. The limiting resistor R_D serves to limit the current, and to make this current approximately constant. To do so, R_D should be much greater than the sample resistance R_{Stot} , which is the total resistance of a sample between the *outside* current contacts shown in the figure above, and also includes the resistance of these current contacts (a contact to a sample may have appreciable resistance, see Appendix B); it is the total resistance encountered by the current as it flows through the sample.

The value of R_{Stot} will be changing in the course of a measurement, as we change the sample temperature. Let us say its maximum value is $R_{Stot} = 500 \text{ k}\Omega$, and the minimum value can be neglected compared to $500 \text{ k}\Omega$, which is indeed the case with VO_2 (assuming contact resistance is not very large). Suppose we decided that we will be satisfied with keeping current constant to 5 %, i.e. we are not worried about up to 5 % possible overall distortion in our measurement. The maximum and minimum rms currents then will be $I_{\max} = V_0/(R_0 + R_D)$ when R_{Stot} is negligible, and $I_{\min} = V_0/(R_0 + R_D + R_{Stot})$ when R_{Stot} is maximum; here V_0 is the rms voltage at the oscillator output. We want to have

$$\frac{\Delta I}{I_{\min}} = \frac{I_{\max} - I_{\min}}{I_{\min}} = 0.05, \quad (2)$$

or, as it is easy to see, simply $R_{Stot}/(R_0 + R_D) = 0.05$, or $R_0 + R_D = R_{Stot}/0.05 = 10 \text{ M}\Omega$. Compared to this value, $R_0 = 50 \Omega$ can be neglected, and essentially our limiting resistor R_D should be set at $10 \text{ M}\Omega$ (which, incidentally, is the maximum resistance which can be set in our decade box). If we use this large limiting resistor, the current will be equal to $I = 2 \text{ V}/10 \text{ M}\Omega = 0.20 \mu\text{A}$.

Let us now discuss the voltage signal. Small current I will result in a fairly small voltage IR_S , especially when R_S also becomes small. Lock-in can measure small voltages, but below about the 1 mV level, the precision will suffer. If we want to increase the current, in the present scheme we will have to sacrifice some additional accuracy in our overall measurement of R_S over the whole range of its variation.

The sample has four contacts: two outside contacts for the current and two inner ones for the voltage. You should read about the four-probe measurement in Appendix B, and a brief description of this principle should be given in your Final Report. Be ready to answer this question: what is the advantage of the four-terminal method over the two-terminal method, and how is this advantage achieved? The voltage from the inner leads V_{meas} (see the figure above) is supplied to the lock-in's measuring input, where it will be amplified (thus the name lock-in *amplifier*), and provided for the readout on the front panel and, through the digitizing interface, to the computer. The resistance between the inner contacts will be found as $R_S = V_{\text{meas}}/I$. Above, when we were discussing the current, we considered R_{Stot} between the outside current contacts, but now we are measuring resistance R_S between the inner voltage contacts. The lock-in will display the RMS value as a DC-level signal.

The maximum voltage V_{meas} in our example (measuring VO_2) is going to be found at room temperature; it will correspond to the maximum value of R_S . Suppose the maximum R_S is equal to $300 \text{ k}\Omega$ (note that we are using here $R_S < R_{Stot}$). Then the maximum

$$V_{\text{meas}} = IR_S = (2 \times 10^{-7}) \times (3 \times 10^5) = 60 \text{ mV}. \quad (3)$$

Such a signal would be quite sufficient, but the trouble is that it will be a lot lower once R_S decreases at high temperatures, where, above the semiconductor-to-metal transition (see below),

we may get a signal as small as about 0.1 mV. The lock-in is going to give us poor resolution in that range. Recall also that, with our proposed numbers, we will have a distortion in resistance values of about 5 % over the whole range, as our current will increase by 5 % when R_S drops by such a large factor.

An alternative way of dealing with changing R_S is to measure the current instead of stabilizing it; it can be best measured by monitoring voltage across the known limiting resistor. Then we do not need such a large limiting resistor: instead of being 10 M Ω it can be, for example, 10 k Ω , and the current will be 1000 times higher, i.e. 0.200 mA instead of 0.200 μ A. But we will have to monitor the current and to use the measured current I_{meas} in a calculation of a sample resistance at every point along the way. This can be done with the use of an additional computer-interfaced voltmeter, and with a modification of PyFERRO.

Do not overdrive the sample. The current should not be increased indefinitely; excessive currents can produce self-heating and even damage the thin-film sample. Do not use currents greater than 1 mA with our VO₂ samples.

So, how is the lock-in measuring resistance? The oscillator output, as we said above, provides sinusoidal voltage of rms value V_0 at a fixed stable frequency which we can set at will. The readout is tuned to read only voltage at this frequency (actually, at a double frequency), all the other frequencies being filtered out. Additionally, it will read only a specified phase of this signal. For our resistive measurement we expect that the specified phase should be zero, as resistor does not produce a phase shift. So, if we are measuring pure resistance, we should be interested in the *in-phase* signal (see Appendix A). The circuit and the phasor diagram of a pure resistor R_S are shown in Figure 4.

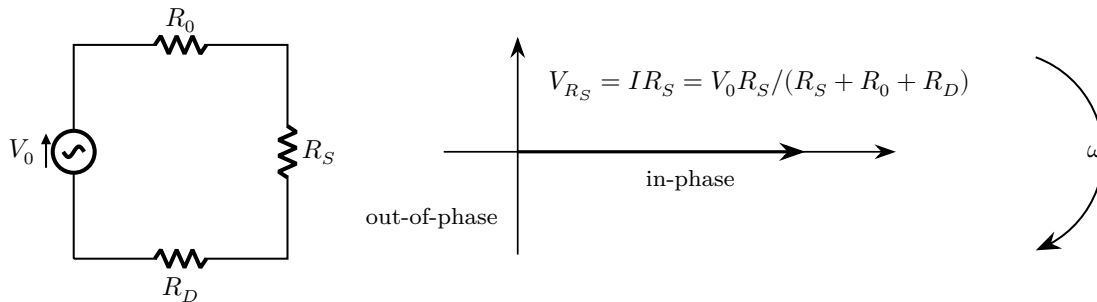


Figure 4: Circuit and phasor diagram for a purely resistive sample.

As always, the phasor should be imagined rotating counterclockwise with ω ; the in-phase direction rotates with it; the projection onto the horizontal axis is the real signal. The current is determined by the source voltage V_0 applied across all of the resistors. We see that in any other direction the signal will be smaller than in this in-phase direction. Thus, if we adjust the phase till we see the *maximum* voltage at the measuring input, this maximum voltage will correspond to the true value of resistance according to $R_S = V_{\text{meas}}/I$. Simultaneously, the 90°-shifted, out-of-phase signal should be zero. It is easier to tune the out-of-phase signal to zero than to find the true maximum in-phase, because when measuring the near-zero signal we can increase the sensitivity. So the usual procedure is to change the phase till one brings out-of-phase signal as close to zero as possible and then to shift the phase by 90° displaying the in-phase signal. Further, all of the above should be independent of frequency $f = \omega/2\pi$. This is what one expects; however, there are complications.

Complications arise because, in addition to the resistance, the circuit contains the so-called *parasitic* or *stray* capacitance C_p . This is the capacitance of shielded cables, etc. In our circuit it

is typically of the order of 10 to 100 pF. This parasitic capacitance can be thought of as being plugged in parallel to the sample resistance R_S , as shown in the circuit diagram above where it is drawn in dashed lines. It will effectively shunt (short out) the resistance, leading to incorrect, decreased values of V_{meas} .

Solution. Recall that lock-in allows you to adjust the oscillator frequency f , which of course implies that reactance $X_C = 1/(2\pi f C_p)$ can be controlled. The solution is then clear: in order to minimize the detrimental effect of C_p we need to decrease f as much as possible: at low f the reactance $1/(2\pi f C_p)$ will be large, and the shunting effect of the parasitic capacitance will be minimized. On the other hand, using extremely low frequency $f \sim 1$ Hz may be difficult, as the lock-in tends to lose its ability to lock the phase at extremely low frequencies. Additionally, the so-called flicker or $1/f$ noise increases at low frequencies. But we do not need to reduce f that much: using $f \approx 50$ Hz is comfortable in terms of this phase-locking and achieves what we need. Indeed, at $f = 50$ Hz, for $C_p = 100$ pF, we find $X_C = 1/(2\pi f C_p) \approx 30$ M Ω , which is more than sufficient to prevent a decrease in V_{meas} due to shunting.

Practical advice. Insert a *known* resistor in place of your sample (preferably, with the value similar to R_S) and empirically find the frequency range at which your measurement agrees with that known resistor value. This is an example of a *baseline* measurement which is done to establish the validity of the measurement technique.

6.2 How to use the lock-in amplifier to measure the capacitance?

The sample is now a capacitor; the material we want to study (BaTiO_3) is a dielectric inside the capacitor. In this measurement we do not need to worry about additionally limiting the current; the current we are dealing with here is sufficiently limited by $R_0 = 50$ Ω . Here we will be using lock-in as a device to measure signals that are *phase shifted* from the phase of V_0 due to the presence of resistors and capacitors in the circuit. The basic circuit diagram for the capacitance measurement is shown in Figure 5.

As in a resistance measurement described above, the signal is generated by the lock-in amplifier's internal reference (oscillator) output; in this measurement we will be setting it at $V_0 = 0.500$ V, $f = 25$ kHz. As before, the built-in internal resistance of the signal generator is $R_0 = 50$ Ω , and there exists parasitic capacitance C_p shown in dashed lines; this capacitance can be thought of as connected in parallel to sample capacitance.

Surprising as it may sound, we will first do the measurement of an *empty* circuit, without the sample. Later we will insert the sample by connecting sample leads 2 and 4 between points 1 and 3, as shown in the diagram above. The purpose of the measurement without the sample is in tuning the lock-in, so as to make it later correctly measure the sample's capacitance.

First, leaving the sample unconnected, we will tune the lock-in so that it is in-phase with the oscillator; it will measure $V_0 = 0.500$ V. We then change the phase by 90° and look at V_{meas} while gradually increasing the sensitivity of the measuring readout. On a sensitive-enough scale we will see some non-zero signal due to parasitic capacitance of our "empty" circuit. We will adjust the phase watching the signal on increasingly more and more sensitive scales till this out-of-phase signal will be made as close as possible to zero. So, in other words, we will follow the same procedure as in measuring resistance, as was described above. See Appendix A for the details of the phasor diagrams in this case.

Once we zeroed the out-of-phase signal, we will insert the sample, so that points 2 and 4 are connected to 1 and 3; once we have done that, the out-of-phase readout from the sample V_{meas} will appear at the measuring input and the magnitude of that out-of-phase signal is going to be proportional to the capacitance of the sample. Indeed, the capacitance of our sample is usually small, of the order of 100 pF–1 nF, and so, with $R_0 = 50$ Ω and at our chosen frequency, we have $\omega R_0 C \ll 1$. In this limit the in-phase component of the signal will be very close to the voltage

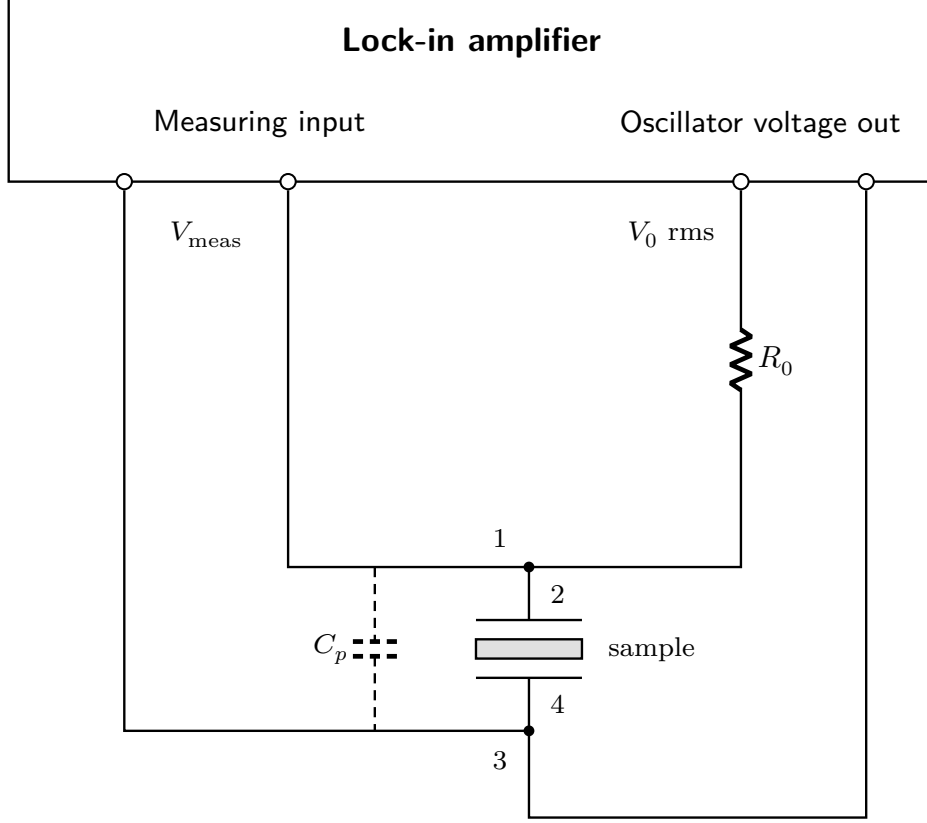


Figure 5: Connecting the lock-in amplifier for a capacitance measurement. $R_0 = 50\,\Omega$ is the built-in internal resistance of the signal generator, and C_p is the parasitic capacitance of the empty circuit, which acts in parallel with the sample.

provided by the signal generator, V_0 , and the out-of-phase component will be close to $V_0 \omega R_0 C$ (see Appendix A for the details). By measuring the out-of-phase component V_{meas} , knowing V_0 , ω and R_0 , we can calculate the capacitance of the sample as

$$C = \frac{\delta}{\omega R_0}, \quad \text{where} \quad \delta = \frac{|V_{\text{meas}}|}{V_0}. \quad (4)$$

6.3 Lock-in amplifier settings

The lab's lock-in amplifier is the Stanford Research Systems SR830 DSP lock-in (Figure 6). Its standard settings for the capacitance measurement are:

Setting	Value
Input	A, AC coupled, shield floating
Sensitivity	1 V
Line notch filters	out
Dynamic reserve	low noise
Oscillator level	0.500 V initially
Oscillator frequency	25 kHz
Reference phase	tuned (see below)
Reference	internal
Oscillator out	SINE OUT to input A (BNC cable)
Time constant	300 ms, 12 dB/octave
Displays	X on channel one, Y on channel two

Setting the SRS SR830. Table 1 gives the same settings in more detail: where to find each one on the front panel, and the equivalent GPIB command. A command followed by a question mark (SENS?) makes the lock-in reply with the current value, which is how PyFERRO reads these settings back and records them in the data-file header. Most of the table is the SR830's own standard set-up, which it returns to if you hold down the Setup key while switching it on; then only the frequency, the level and the time constant need changing. PyFERRO can also send the whole table (Appendix C, Lock-in tab).

Oscillator resistance and signal size. The SR830's sine output has an output resistance of $R_0 = 50\ \Omega$ (marked under the SINE OUT connector). The lock-in's least sensitive range is 1 V, and the in-phase signal is close to the oscillator level V_0 , so V_0 must stay below 1 V; at $V_0 = 0.500\text{ V}$ it is half of full scale. The out-of-phase signal is then small: for 100 pF at 25 kHz, $V_{\text{meas}} = V_0 \omega R_0 C \approx 0.4\text{ mV}$. Read it with the channel-two Expand (below). A known resistor R_s in series with SINE OUT raises the signal in proportion; then use $R_0 = 50\ \Omega + R_s$, and check that $\omega R_0 C$ stays well below 1 or use the exact formula of Appendix A.

The SR830 is a dual-phase lock-in: it shows the in-phase signal as X and the out-of-phase signal as Y at the same time. Tune the phase with the sample disconnected until Y is as close to zero as possible (Auto Phase does this), then connect the sample and read V_{meas} from Y .

Before measuring the sample, we do a practice run by measuring the capacitance of a piece of BNC cable. See Appendix A for the relevant phasor diagrams.



Figure 6: The SRS SR830 DSP lock-in amplifier. Photo: Stanford Research Systems.

Setting	Value	How to set it on the front panel	GPIB
Reference	internal	Press Source, in the REFERENCE section at the right, until INTERNAL is lit.	FMOD 1
Oscillator frequency	25.000 kHz	Press Freq and turn the knob; the REFERENCE display shows the frequency.	FREQ 25000
Oscillator level	0.500 V	Press Ampl and turn the knob.	SLVL 0.5
Signal input	A	Cable to the A/I input. Press Input, in the SIGNAL INPUT section at the lower left, until A is lit.	ISRC 0
Coupling	AC	Press Couple until AC is lit.	ICPL 0
Grounding	float	Press Ground until FLOAT is lit.	IGND 0
Line notch filters	out	Press Notch, in the FILTERS section, until neither LINE nor 2xLINE is lit.	ILIN 0
Sensitivity	1 V	The arrow keys in the SENSITIVITY section; the up arrow moves to a larger full scale (less sensitive).	SENS 26
Dynamic reserve	low noise	Press Reserve, in the RESERVE section, until LOW NOISE is lit.	RMOD 2
Time constant	300 ms	The arrow keys in the TIME CONSTANT section.	OFLT 9
Filter slope	12 dB/oct	Press Slope/Oct, in the TIME CONSTANT section, until 12 dB is lit.	OFSL 1
Offset and expand	off	Under each display: OFFSET On/Off until the offset is off, and Expand until neither x10 nor x100 is lit.	OEXP 1,0,0 (and 2, 3)
Displays	X and Y	Press Display under each display until X (channel one) and Y (channel two) are lit.	DDEF 1,0,0, DDEF 2,0,0
Reference phase	Y = 0 with the sample disconnected	Press Phase in the AUTO section (Auto Phase); repeat once the outputs have settled. For fine tuning, press Phase in the REFERENCE section and turn the knob.	APHS

Table 1: SRS SR830 settings for the capacitance measurement. The GPIB commands are from chapter 5 of the SR830 manual [1]; **SENS 26** and **OFLT 9** are positions in its tables of sensitivities (2 nV to 1 V) and time constants (10 μ s to 30 ks).

First, check the signal generated by the lock-in amplifier’s internal oscillator. Connect one end of a BNC cable to SINE OUT and the other end to the A/I input, and set the lock-in as in Table 1. Press Phase in the REFERENCE section and turn the knob to set the reference phase to 0.00°. Channel one shows the in-phase component X and channel two the out-of-phase component Y. If most of the signal appears in X, its reading should be close to the oscillator setting of 0.500 V.

We are interested in the out-of-phase component Y. Ideally, Y should be zero after the phase has been correctly adjusted. However, the electronics inside the lock-in amplifier can introduce a small phase error, so Y may initially have a small but nonzero value.

The lock-in amplifier provides an automatic way to correct this phase error. Press Phase in the AUTO section to run Auto Phase. Wait a few seconds for the readings to settle. The in-phase component X should then be maximized, while the out-of-phase component Y should be close to zero.

Upon closer inspection, you may find that Auto Phase has left a small out-of-phase signal, too small to read on the 1 V range. Keep the sensitivity at 1 V and press Expand under channel two once (x10) or twice (x100). The expand does not change the value shown for Y: it adds digits to it, so a signal of a fraction of a millivolt can be read. Press Phase in the REFERENCE section and turn the knob until Y is as close to zero as possible. If channel two’s display shows OVLD, reduce the expand. Do not reduce the 1 V sensitivity setting while the oscillator is connected directly to the input: X would then be beyond full scale.

Now that the out-of-phase signal has been brought as close to zero as possible, keep the sensitivity at 1 V. Insert a BNC “T” connector into the line and connect another BNC cable. The added cable will produce a nonzero out-of-phase signal V_{meas} , which you can read from Y. Calculate its capacitance using $C = |V_{\text{meas}}|/(V_0\omega R_0)$, where $V_0 = 0.500 \text{ V}$, $\omega = 1.57 \times 10^5 \text{ s}^{-1}$, and $R_0 = 50 \Omega$. Measure the length of the cable, then repeat the measurement with a longer cable. Turn the expand off again when you are done: PyFERRO notes an expand in the data file.

Question. In round numbers, what is the capacitance of one foot of BNC cable?

7 The rest of the experimental apparatus in detail

7.1 Sample chamber

A photograph of the sample chamber is shown in Figure 7. There are four wires reaching the inside of the chamber, “red” and “black” pairs. The sample will be either a film of VO_2 with four resistive leads, or a capacitor with BaTiO_3 requiring only two leads attached to its plates. If two students do the measurement, it is recommended that each one of them mounts their own sample and evaluates the data independently. However, VO_2 samples are fragile, and a staff member should be asked to mount them in the chamber, at least initially.

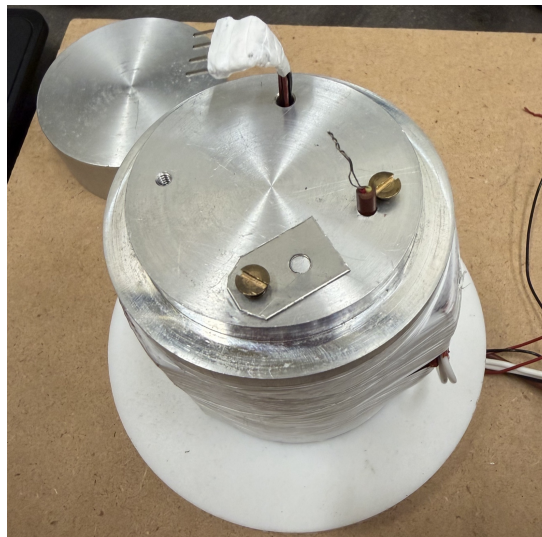


Figure 7: The sample chamber with its metal cup cover removed, showing the pins for connecting the PCB and the cable to the PID temperature controller.

In normal use, the sample chamber is covered by an aluminum cup (see the upper left corner of the photograph). The whole assembly is covered by a stainless steel cover (the dewar) for thermal insulation. If you wish to cool the sample faster, remove the dewar. You may also remove the aluminum cup, but in that case stop taking data, since the temperature readout will not be reliable.

7.2 Samples and sample preparation

Samples of VO_2 in a form of thin films (about 100 nm thick) have been deposited on Si/SiO₂ substrates in another lab on campus. Contacts to them will be attached by staff members. These samples are rather fragile, especially those contacts, and should be handled with great care. Ask a staff member to plug the sample into the holder, at least initially.

The BaTiO₃ sample is from a commercial distributor. It comes in ceramic form with many different, randomly oriented microcrystals. Feel free to split off your sample from larger chunks of the ceramic blocks. Pick a flat and thin piece of ceramic BaTiO₃, and shape it to a thin plate by lightly sanding on the sides. Keep in mind that capacitance $C = \epsilon A/d$: the thinner the better, and the larger the area the better. Also remember that you will need to estimate A and d numerically, and so the more regular the sample shape, the easier and more accurate this estimate will be.

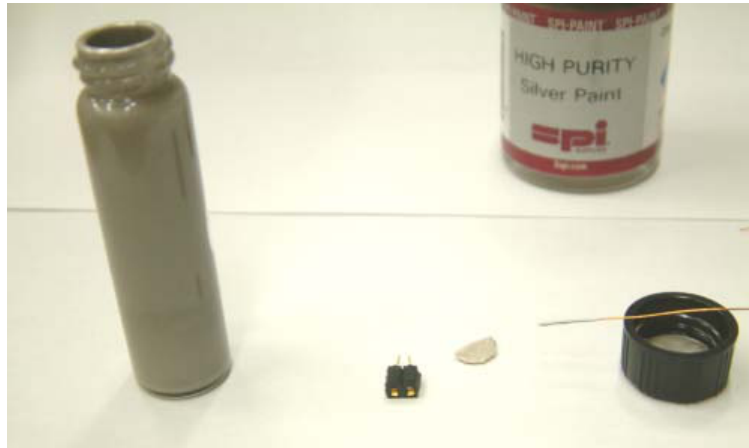


Figure 8: Materials for mounting the BaTiO₃ sample: silver paint, a two-pin connector, a ceramic chip, and the thin wire used to apply the paint.

The procedure for preparing and mounting the BaTiO₃ sample is as follows.

1. Sand the BaTiO₃ sample to a uniform size and measure its dimensions.
2. Apply silver paint to its two broad faces to form the capacitor plates, then sand the unpainted sides so that the two plates do not touch.
3. Check for a short circuit with an ohmmeter.
4. Apply a thin layer of silver paste to the PCB pad connected to one terminal of the through-hole connector, and place one painted face of the sample onto the paste.
5. Solder a wire to the pad connected to the other terminal and connect the wire to the painted upper face of the sample. The two painted faces are thus connected to separate PCB pads and separate connector terminals.
6. Check for a short circuit again through the two connector terminals. Be careful not to let the silver paint or paste bridge the two faces: you are making a capacitor.
7. Install the sample in the chamber, the entire PCB is secured in place by the metal clip. The two black two-pin connectors on the PCB plug into the first two connector pins in the sample chamber, providing the electrical connections to the two painted faces of the sample. You may have to check which lead goes where with a voltmeter.

The paint dries very quickly. Do not leave the silver-paint container open for an extended period of time; open it for as short an interval as possible, just long enough to take the paint you need.

A stereo microscope is provided for mounting the sample, although with good eyes and eye-hand coordination it can be done without the microscope.

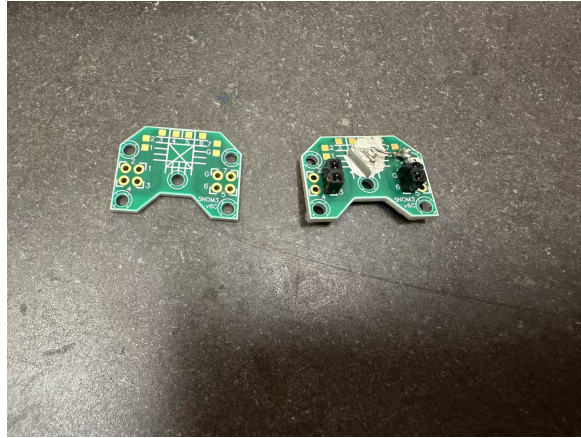


Figure 9: An unloaded sample-mounting PCB (left) and a PCB with a mounted BaTiO₃ sample and two two-pin connectors (right).

Return your sample. The sample of BaTiO₃ (but not of VO₂) you make for the studies should be returned with your final lab report. Simply attach the sample with Scotch tape to the first page of your report.



Figure 10: Wiring between the sample chamber and the four-connector BNC panel.

Connect the experimental apparatus: one BNC cable from the signal source to input A of the lock-in, and another cable from the lock-in to the BNC leading to the sample. Check that

there is no sample in the sample chamber. (For later use, when you put in the sample, you need to terminate the other BNC with a short circuit to ground.)

The Y reading (the out-of-phase component) that you see now is due to all of the capacitances leading to the sample, including the capacitance between the Teflon-insulated black wire and ground. Fine-tune the phase so that this background Y reading is as close to zero as possible. Keep the sensitivity at 1 V and insert the sample. The additional Y signal is due primarily to the sample; use its measured voltage V_{meas} in $C = (V_{\text{meas}}/V_0)/(\omega R_0)$ to calculate the sample's capacitance. If Y is too small to read clearly, expand the Y display as described in Section 6.3, and turn the expand off again before recording.

For all of this to work, $V_{\text{meas}}/V_0 = \omega R_0 C$ should be much less than 1. In Appendix A we derived the exact formula for calculating the capacitance, and noted that the difference between the exact and approximate formulae is about 1 % when $\delta = 0.1$. Therefore, for an initial oscillator setting of $V_0 = 0.500$ V, the systematic error caused by using the approximate formula is about 1 % or less if V_{meas} is smaller than 50 mV. If the measured oscillator voltage differs from this setting, use the measured value of V_0 and adjust the voltage limit accordingly. Here $R_0 = 50 \Omega$ and, at 25 kHz, $\omega = 2\pi f \simeq 1.57 \times 10^5 \text{ s}^{-1}$. If you made a very thin sample and its capacitance is large, you may go outside this limit. In that case, use the exact formula or reduce the frequency to 10 kHz or 5 kHz and repeat the phase tuning. If you change the frequency, use its new value of ω when calculating the capacitance.

At this point you are ready to measure the temperature dependence of the capacitance. Keep in mind that we expect the capacitance to increase as we warm the sample, and select a range where the current signal is less than one third of the full range. This way you do not have to change sensitivity when the most interesting part of the measurement is happening.

7.3 Temperature control

The front panel of the Omega CN08D3-SR-2U-CM-AC temperature controller [4] is shown in Figure 11. The controller is connected to a nearby solid-state relay. Its low-power control output turns the relay on and off; the relay, in turn, switches the separate power circuit that drives the heater wrapped around the sample chamber (see the photograph of the sample chamber above). In PID operation, changing the fraction of time that the relay is on adjusts the average power delivered to the heater. The input for the temperature controller is a thermocouple that is placed in the aluminum block of the sample chamber. A thermocouple is a pair of wires made of two different metals, joined at one end. When the joined end is at a different temperature from the free ends, a small voltage appears between the free ends (the Seebeck effect), typically tens of microvolts per degree Celsius.

The simplest form of temperature control is “on/off” control. You specify a *setpoint*, the highest temperature needed for the measurement. When the measured temperature is below it, the controller turns the heater fully on; when the setpoint is reached, it turns the heater off. Heat remaining in the chamber may continue to raise its temperature after the heater turns off. This is called “overshoot”. As the temperature falls below the setpoint, heating starts again. Because different parts of the chamber take time to warm and cool, the temperature may continue to fall before rising again and may oscillate around the setpoint. “PID” control adjusts the heating output to reduce these oscillations; see Ref. [4] for an introduction to PID control. Compare the temperature response under on/off and PID control.

For a phase-transition measurement, the heating rate also matters. In the CN08D3's Slope SV (SLOP) mode, you specify a target temperature and a rate at which the setpoint rises toward it. PID control adjusts the heater output, reducing it when necessary to limit the rate of heating. This helps the sample temperature follow the measured chamber temperature and is important for obtaining reliable data. Try different heating rates, examine the data, and choose a suitable rate based on their quality. The CN08D3 also has a Program SV (PROG) mode for scheduling multiple heating and holding steps, but it is optional for this experiment.



Figure 11: Front panel of the OMEGA CN08D3-SR-2U-CM-AC temperature controller, showing the measured temperature and setpoint.

Question. Be ready to answer this question: what is the role of the “I” in PID?

Front-panel operation. The controller is operated using the gear-symbol key, the curved-arrow key, and the down and up arrow keys. From the normal PV/SV display, a brief press of the gear-symbol key enters Regulation Mode; holding it for more than three seconds enters Initial Setting Mode. The curved-arrow key moves through the menu items, and the arrow keys change a value or option. Press the gear-symbol key to save a setting. The F1 and F2 keys have user-defined functions and are not needed here. On the normal display, PV is the measured temperature and SV is the setpoint.

To set the temperature unit to Celsius, hold the gear-symbol key for more than three seconds to enter Initial Setting Mode. Use the curved-arrow key to find $tPUN$, then use the up or down arrow to select $^{\circ}C$. Press the gear-symbol key to save the setting and return to the normal display. Confirm that the display shows $^{\circ}C$, and check that SV is set to the desired temperature.

To select slope heating, hold the gear-symbol key for more than three seconds to enter Initial Setting Mode. Press the curved-arrow key until CRTS appears, use the up or down arrow to change CONS to SLOP, and press the gear-symbol key to save and return to the normal display. The slope-rate unit is set separately in Initial Setting Mode using EVEC: its third digit from the left is 0 for $^{\circ}C\ min^{-1}$ or 1 for $^{\circ}C\ s^{-1}$. Change only this digit if the rate unit needs to be changed. To set the heating rate, briefly press the gear-symbol key to enter Regulation Mode, use the curved-arrow key to find SVSL, set the desired rate with the up or down arrow, and press the gear-symbol key to save, as shown in Figure 12.

To begin heating, press the curved-arrow key from the normal display until R-S appears. Select RUN with the arrow keys and press the gear-symbol key to save. Monitor PV and the measurement data during heating. When the target temperature is reached and cooling data collection is to begin, use the curved-arrow key to find R-S again, select STOP, and press the



Figure 12: The OMEGA CN08D3-SR-2U-CM-AC front-panel display for setting the SV slope (heating rate).

gear-symbol key to save. Confirm that the heating output is off, then continue recording as the sample cools. In **SLOP** mode, reaching the target temperature does not automatically stop the controller.

Hot surfaces and overheating. Most of the capacitance measurements are made at temperatures over 100 °C. Be careful: the sample chamber will be hot. Do not exceed 160 °C, because the sample chamber may be damaged. The sample chamber can be seriously overheated if left unattended with the heater on. Always switch off the temperature controller when the measurement is done. When using **SLOP** mode, reaching the target temperature does not turn off the heater. Set **R-S** to **STOP** and confirm that the heating output is off before recording the cooling data.

7.4 Temperature measurement

The temperature is measured by a thermocouple installed in the aluminum block of the sample chamber and connected to the temperature input of the Omega CN08D3 controller. A thermocouple detects temperature through the small voltage between its two wires, which grows with the temperature difference between the junction in the block and the controller's terminals.

The CN08D3 measures the thermocouple voltage, corrects for the temperature of its own terminals, and converts the result to the process temperature (PV). Its input-type setting must match the thermocouple type. PyFERRO reads the PV and setpoint (SV) from the controller over RS-485 and records them with the lock-in measurements. The controller uses the PV as feedback to regulate the heater. Because the probe measures the aluminum block, its reading may differ slightly from the temperature at the sample.

The PyFERRO data file contains (1) the time dependence of the temperature and (2) the temperature dependence of the lock-in output signals *X* and *Y*. Copy the data file to a thumb drive for analysis.

Before measuring the sample, use PyFERRO to study the temperature control process. On the Instruments tab, turn off Simulation mode and check the connections to the lock-in amplifier and the Omega CN08D3 temperature controller. On the Run tab, select the CN08D3 controller as the temperature source. Start at room temperature or below 50 °C. On the CN08D3, select constant setpoint mode (CONS) and set the setpoint to 70 °C. Set the controller's CONTROL mode to ON-OFF. Start monitoring and recording in PyFERRO, then set the controller to RUN. Continue recording as the temperature rises and oscillations develop. Then change the CN08D3's CONTROL mode from ON-OFF to PID without stopping the recording. Compare the temperature-versus-time data for the two modes and discuss your observations in your lab report.

8 Suggested protocol for a measurement

1. Tune the lock-in amplifier and place the VO₂ or BaTiO₃ sample in the sample chamber. Close the chamber with its fitted cover. If additional thermal insulation is needed, place the larger stainless-steel dewar over the covered chamber.
2. Start the PyFERRO routine (see Appendix C).
3. Set the temperature-controller setpoint on the front panel to 90.0 °C for a VO₂ measurement or 150.0 °C for the capacitance measurement of BaTiO₃.
4. Program the ramp rate. Change the ramp rate and examine the data: you want to make sure that your ramp rate is slow enough for the data to be reliable. Discuss with the laboratory staff what this means.
5. Set R-S to RUN to begin heating, the green light on the Solid-State Relay should blink.
6. Watch the signal as the sample crosses the phase transition. In VO₂ this will take place over a fairly wide temperature range; in BaTiO₃, in a narrow temperature range.
7. Once the setpoint set R-S to STOP on the temperature controller and confirm that the heating output is off. You will take more data while the sample is cooling.
8. The sample will slowly cool while the data is recorded. In the BaTiO₃ measurement, nothing interesting happens below 80 °C, so when the temperature reaches 80 °C you can take off the dewar to speed up the cooling. When measuring a VO₂ sample, you want to take the data in a slow-enough ramp mode in the whole temperature range.

9 Advice for evaluating the capacitance data

From the graph of the measured capacitance as a function of temperature you will find the T_C for this material.

Estimate the area of the capacitor and the distance between the plates and calculate the relative dielectric constant k at room temperature. To find k , you will need to calculate the capacitance C_0 of an empty (i.e. filled with vacuum or air) capacitor with the same geometry as your sample, and then take your measured C and divide it by that calculated C_0 : $k = C/C_0$.

Give an estimate of the error in k due to the necessarily inaccurate measurement of the area and the distance between the painted plates. This error is expected to be large. Compare the value of k to the known dielectric constant of this ceramic material.

Our main interest, however, is to use the data from our ceramic sample to access some of the intrinsic, fundamental properties of BaTiO₃. This means that your values of k will need an important correction, which we discuss here. If the dielectric constant is very large, the capacitance of the sample should be very large too. However, in this measurement, the

capacitance will not increase near the transition that much. A possible explanation is that in a ceramic sample the dielectric polarization is interrupted at the grain boundaries. The capacitance between the grains provides coupling, but that inter-grain capacitance does not change with temperature, and certainly does not diverge at T_C .

You can model this behavior by taking the sample as *two capacitors in series*: C_b (capacitance of the grain boundaries) and C_{BTO} (capacitance of the bulk grains of BaTiO_3). For the effective (measured) capacitance in this case we have

$$C = [C_b^{-1} + C_{BTO}^{-1}]^{-1}. \quad (5)$$

Here C_b is independent of the temperature and C_{BTO} corresponds to an ideal BaTiO_3 uniformly-filled capacitor. In this model you can assume that at the phase transition point, when C_{BTO} is very large, the C_{BTO}^{-1} term can be neglected, and all you measure is C_b , i.e. $C \approx C_b$. Thus, from your measurement of the peak capacitance at T_C , you will know C_b . For all the other temperatures you can calculate C_{BTO} as

$$C_{BTO} = \frac{C_b C}{C_b - C} = \frac{C}{1 - C/C_b}, \quad (6)$$

where C is the measured capacitance and C_b is the maximum (peak) capacitance at T_C . Therefore, you should use this corrected C_{BTO} and not the measured capacitance C directly in finding k values: $k = C_{BTO}/C_0$.

Once you obtain corrected k versus T data, refer to Appendix E (which you should be reading continually while doing this experiment; it contains most of the physics), analyze Curie–Weiss dependence of k in the vicinity of the transition, and find Curie–Weiss constant C_{CW} . From it, calculate the microscopic dipole moment and the displacement of Ti ion in the unit cell based on your data. Compare it to the displacement you will find in the literature (you will have to find this value, which is known from X-ray measurements of BaTiO_3). Also calculate the same quantities using the intrinsic (single-crystalline) value of C_{CW} from the literature.

A Elements of AC circuit theory relevant to our measurements using the lock-in amplifier

We are dealing with sinusoidal AC voltages and currents, such as $V = V_0 \cos(\omega t + \gamma)$, where V_0 is the amplitude, ω is the angular frequency, and $V(t = 0) = V_0 \cos \gamma$ is the initial condition, i.e. a signal at $t = 0$; γ is called the *phase angle*, or the *phase*. You will recall that on a *phasor diagram* which you studied in your first physics course, all signal amplitudes are represented as rotating vectors; ω is the angular velocity of a rotating phasor related to regular frequency by $\omega = 2\pi f$; projections of these vectors onto a chosen axis are real instantaneous signals. A *cosine* function used above corresponds to a projection onto a horizontal axis; if we used a *sine* function, it would represent a projection onto a vertical axis. In Figure 13 we show V_0 at the moment of time t , when it makes an angle $\omega t + \gamma$ with the horizontal; we see V_0 's projection onto the horizontal direction, $V_0 \cos(\omega t + \gamma)$, which is the real signal at the moment t .

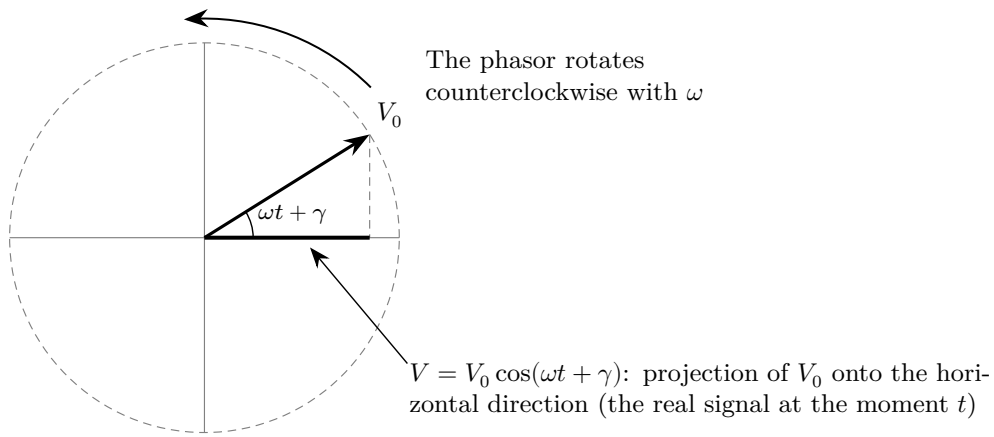


Figure 13: A rotating phasor and its projection onto the horizontal axis.

We will recognize that a phase angle also specifies a *direction* on a phasor diagram; there is an infinite number of possible radial directions in a circle.

It is convenient to measure all phases with respect to the phase of the source, V_0 . A zero phase angle then corresponds to the direction of V_0 ; we will call it in what follows the *in-phase* direction. A direction perpendicular to the in-phase is called *out-of-phase*. It is clear that projections of V_0 onto the in-phase and out-of-phase directions are V_0 and zero respectively.

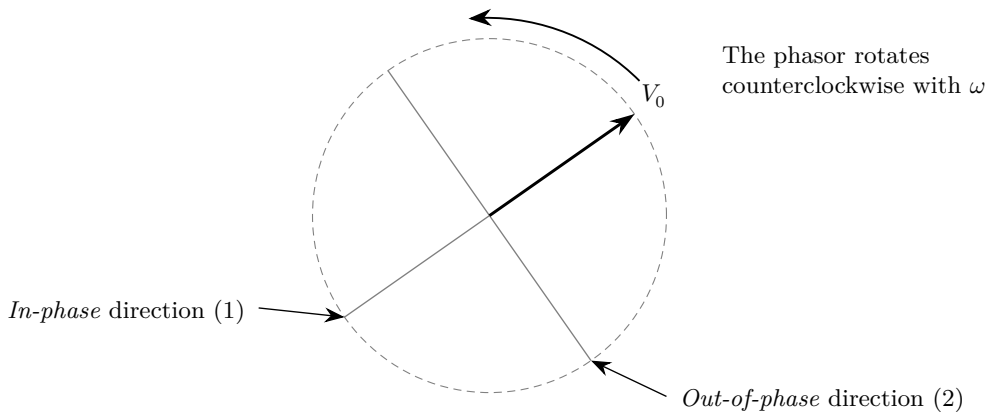


Figure 14: The in-phase and out-of-phase reference directions, defined relative to the source voltage V_0 .

Deciding to use these two phases (or these two phase directions), we are effectively setting our initial γ to zero. But we will be eventually measuring average rms values which do not depend on the initial phase; at the same time, all *relative* phases will be preserved in this new system of phase referencing.

Measuring a specified phase of a signal (or “looking” at a phase, as we say in electronics) is equivalent to looking at a signal projection onto a chosen direction in a phasor diagram. While we chose two such directions, in-phase and out-of-phase with respect to the source voltage V_0 , the lock-in is capable of looking at any phase direction. This is one of its advantages over a normal voltmeter.

With this, let us now consider a phasor of an RC circuit. Figure 15 shows a phasor of a series connection of R and C , with AC source supplying voltage V_0 ; it is a vector sum of voltages V_R and V_C . For completeness we also show the current I , which is established in the circuit; it is shown pointing in the direction of V_R . The current is shown only once because it is the same for all elements in a series connection.

Note that current is said to be *in-phase with the voltage across the resistor*, V_R . We should be careful not to mix our in-phase direction along V_0 as it was defined above with the in-phase direction for the current; these are two different phase directions, as can be seen on the diagram.

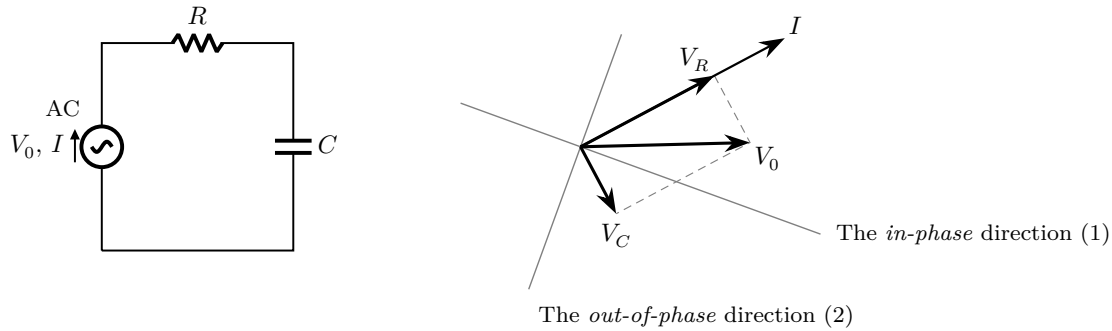


Figure 15: Series RC circuit and its phasor diagram.

It is clear that projections of V_R and V_C onto both in-phase and out-of-phase directions are non-zero, while the algebraic sum of projections of vectors V_R and V_C onto the in-phase direction is V_0 and to the out-of-phase direction is zero. Projections of V_R and V_C constitute signals measured with the lock-in amplifier across R or C respectively, either in-phase or out-of-phase.

Measuring signals in a particular phase as described above is called “phase-sensitive” detection. We are not discussing here how it is practically implemented in the inside electronics of a lock-in amplifier, but only its meaning in terms of a phasor diagram.

Measuring resistance is discussed in the main text in Section 6; we need to watch for the detrimental effect of a parasitic capacitance, and the practical solution is to measure at low-enough frequency.

A.1 Measuring capacitance

Here is an illustration of the lock-in tuning process in terms of phasor diagrams.

Initial tuning. The source voltage V_0 is directly measured, $V = V_0$. The phase is adjusted so that there is no out-of-phase component (Figure 16).

The capacitor is added to the circuit. The voltage drop on the resistor $V_1 = V_R$ and the voltage drop on the capacitor $V_2 = V_C$, via vector addition, add up to the source voltage V_0 .

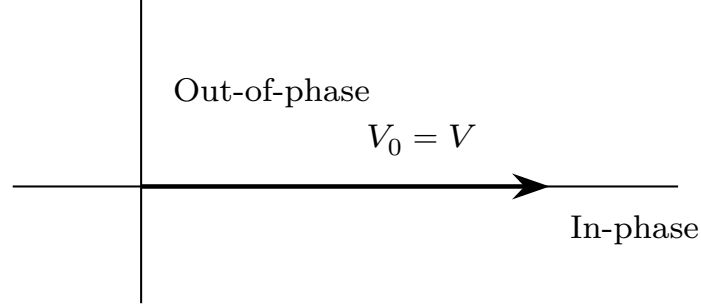


Figure 16: Initial tuning: the measured voltage lies entirely along the in-phase direction.

There is a 90° phase difference between V_R and V_C . The phasor diagram looks as shown in Figure 17. The phase angle α is given by

$$\tan \alpha = \frac{V_R}{V_C} = \omega RC. \quad (7)$$

For $\omega = 2\pi f = 2\pi \times 25 \text{ kHz}$, $R_0 = 50 \Omega$ and $C = 100 \text{ pF}$ we get $\tan \alpha = 0.00079$. For such small angles (in radians!) we can safely take $\tan \alpha = \alpha = \omega RC$.

For $V_0 = 0.500 \text{ V}$ the measured out-of-phase voltage is $V_2 = V_C \sin \alpha$. For small α , $V_C \approx V_0$ and $V_2 = -V_0 \omega R_0 C = -0.39 \text{ mV}$.

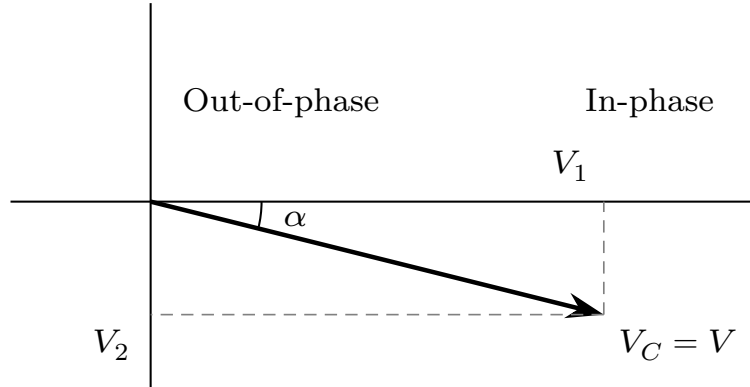


Figure 17: With the capacitor added, the measured voltage acquires a small negative out-of-phase component.

The parasitic capacitance, retuned. Connecting the measurement circuitry is like adding a capacitor to the circuit. Here $C = C_p$ represents the sum of all parasitic capacitances. Now we tune the lock-in so that the out-of-phase component is zero again. The angle α is still small, so $\tan \alpha = \alpha = \omega R_0 C_p$ (Figure 18).

The sample is added. Finally we add the sample. The capacitance gets larger: $C' = C_p + C$, where C is the capacitance of the sample. Since the larger capacitance results in smaller impedance, the current in the circuit increases, and the voltage drop on the resistor increases. The phasor diagram now looks as in Figure 19, where $\tan \beta = \omega R_0 C'$.

The measured out-of-phase signal is $V_2 = V_{C'} \sin \varphi$, where $\varphi = \alpha - \beta$. As long as all of these

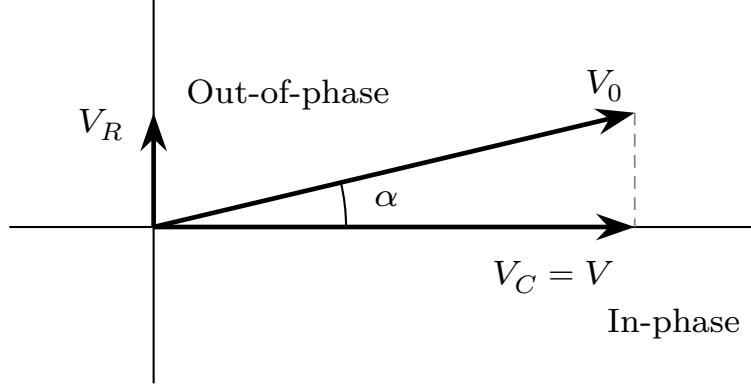


Figure 18: After retuning, the measured signal is again purely in-phase; the source phasor V_0 now sits at the small angle α .

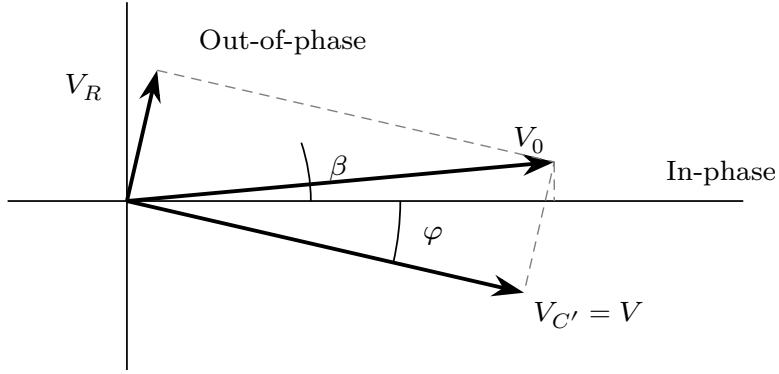


Figure 19: With the sample inserted, the measured out-of-phase signal is $V_2 = V_{C'} \sin \varphi$ with $\varphi = \alpha - \beta$.

angles are small, we get $V_2 = V_0 \omega R_0 C$, where C is the capacitance of the sample. Therefore

$$C = \frac{\delta}{\omega R_0}, \quad (8)$$

where $\delta = -V_2/V_0$ (δ is a positive quantity and, of course, so is the capacitance). Here $V_2 = V_{\text{meas}}$.

The capacitance of the typical sample is in the 100 pF range. For 100 pF, at 25 kHz and $R_0 = 50 \Omega$, the quantity $\omega R_0 C$ is much less than 1. Therefore, in this limit,

$$V_1/V_0 \sim 1, \quad (9)$$

$$V_2/V_0 \sim -\omega R_0 C, \quad (10)$$

$$C = -\frac{V_2/V_0}{\omega R_0} = \frac{\delta}{\omega R_0}. \quad (11)$$

A.2 Alternative derivation using complex numbers

The same result can be obtained using complex numbers instead of phasors. Complex numbers and symbols for treatment of AC circuits were introduced in the “Primer to Linear Circuits and Transmission Lines”, posted with the Course Documents on Blackboard. The impedance of the capacitor is $Z = 1/(j\omega C)$. The current in the RC circuit shown in Figure 20 is $I = V_0/[R + 1/(j\omega C)]$, and the voltage on the capacitor is

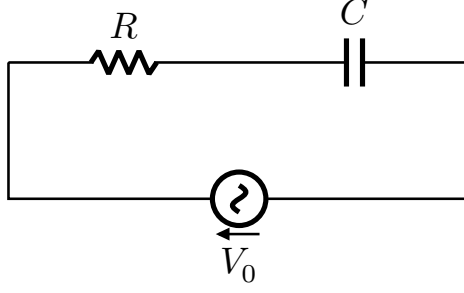


Figure 20: The series RC circuit used in the complex-number derivation.

$$V = \frac{V_0 / j\omega C}{R + 1/j\omega C} = \frac{V_0}{1 + j\omega RC} = \frac{V_0}{1 + (\omega RC)^2} - \frac{j V_0 \omega RC}{1 + (\omega RC)^2}. \quad (12)$$

The real and imaginary parts are

$$V_1 = \frac{V_0}{1 + (\omega RC)^2} \quad \text{and} \quad V_2 = -\frac{V_0 \omega RC}{1 + (\omega RC)^2}. \quad (13)$$

With the proper setting of the lock-in amplifier, described in the main text, we can measure V_2 . Since we know V_0 , ω and R , we can calculate the capacitance of the sample:

$$C = \frac{1}{\omega R} \frac{1 - \sqrt{1 - 4\delta^2}}{2\delta}, \quad (14)$$

where $\delta = -V_2/V_0$ (since the measured voltage is negative, δ is a positive number).

So, to summarize, when $\omega R_0 C$ is small ($\ll 1$) we can use formula $C = \delta/(\omega R)$, but for a very large and thin sample, capacitance will not be small, and one has to use the complete expression (14) derived above. A simple calculation shows that the difference between the exact and the approximate formula is about 1 % when $\delta = 0.1$ and about 5 % when $\delta = 0.2$.

B The four-probe method: advantages and limitations

A simple resistance measurement is performed either by measuring both the current through the resistor and the voltage across the resistor (an I and V measurement), as in Figure 21,

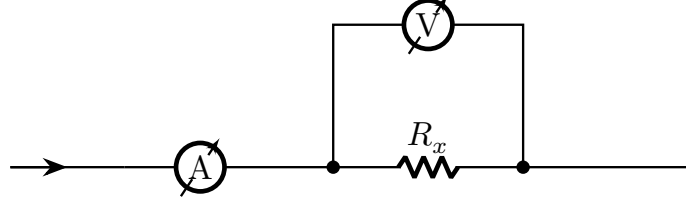


Figure 21: Two-terminal resistance measurement.

or by measuring a current, supplied by a battery, which flows through an adjustable resistor and the unknown resistor connected in series, with the scale being set as to read the unknown resistance (an *ohmmeter* measurement, which we do not show here). Either way, it is a *two-terminal* measurement: the unknown resistor is connected to the rest of the circuit at two points, or two terminals, shown in Figure 21 as black dots. In this measurement it is implicitly assumed that *contact resistance* and resistance of the wires is negligible compared to the resistance R_x . In many cases, especially when dealing with commercially-made resistors, this is an excellent assumption.

There are however many experimental situations in which this measurement is complicated by the presence of finite, non-negligible *contact resistance*, and sometimes also by the presence of *lead resistance* (resistance of electrical connections, wires). Consider for example a measurement performed on a sample shaped as a uniform-width strip (for example, it may be a thin film sample). When performing a two-terminal measurement, we place two metal contacts at the ends of a film strip and attach to them current and voltage leads, as shown in Figure 22.

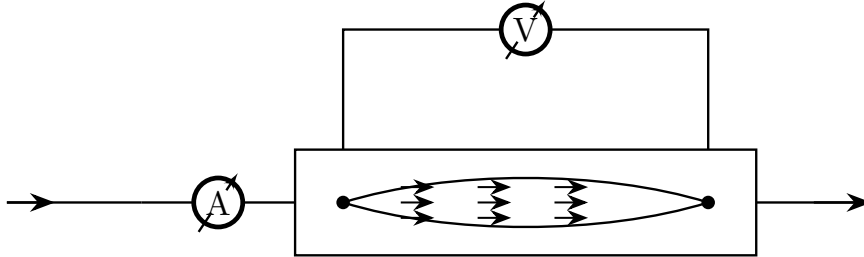


Figure 22: Two-terminal measurement on a sample shaped as a strip, with two contacts.

There are at least two problems with this arrangement:

1. If there is appreciable contact resistance between the film and the contacts (black dots), the equivalent electrical circuit describing this measurement will look like the one depicted in Figure 23, where the resistors R_{cont} indicate the presence of said contact resistance. If R_{cont} is not negligibly small compared to R_x , the voltage drops across them, $2IR_{\text{cont}}$, will be included in the voltmeter measurement. Note that here the full sample current I flows through these R_{cont} .
2. Current flow lines in Figure 23 may be non-uniform around the contacts. This makes the length of a sample ill-determined, and thus will introduce an error in the calculation of

resistivity, even if contact resistances were negligible.

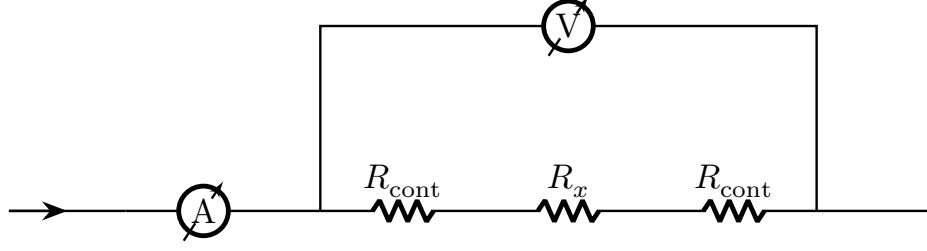


Figure 23: Equivalent circuit showing the presence of contact resistance in a two-terminal measurement.

We note that the same un-avoidable problems exist in an ohmmeter two-terminal measurement.

These problems with the two-terminal connection can be overcome by making a *four-terminal measurement* described below. This measurement is known under different names such as the *four-probe method*, a *four-terminal sensing*, a *four-probe geometry*, or as *Kelvin sensing* (the latter name is rarely used today; it acknowledges that a similar method was first described by Lord Kelvin a very long time ago). It is universally used in solid state physics to perform accurate measurements of true material parameters, primarily of resistivity ρ , as well as in other situations when contact resistance needs to be excluded from the measurement.

The idea is simple: make current and voltage contacts separately. The current is measured, and so the results will not depend of current contact resistance. The two inner voltage leads detect *potential difference* between them. The current practically does not flow into these inner leads, as the voltmeter resistance is very high. So again, their contact resistances do not enter into the picture, except for the case when these resistances are becoming very large, namely, comparable with the input resistance (impedance) of a voltmeter.

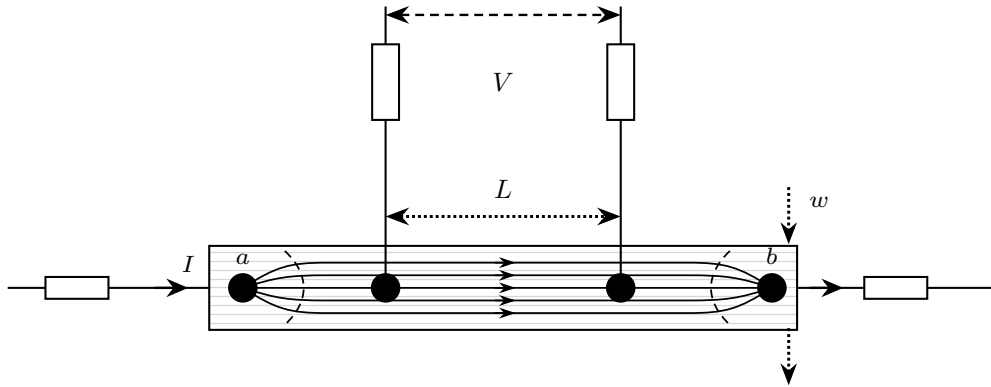


Figure 24: The basic four-probe method.

If the current is I and the potential difference between the inner leads is V , then the resistance of the sample which lies between the voltage leads is simply $R_S = V/I$. Now we can find the resistivity ρ from the formula

$$\frac{V}{I} = R_S = \frac{\rho L}{w\delta} \quad \text{or} \quad \rho = R_S \frac{w\delta}{L}, \quad (15)$$

where w is the sample width, δ is its thickness, and L is the distance between the voltage leads. Here we assume that we know or can measure sample geometry (uniform width w , uniform thickness δ , distance between voltage leads L).

As we said, this method was designed to exclude from the measurement contact and wire resistances, and it always achieves this goal once we separate the two sets of contacts. However, the question of *how accurately* one can find the absolute value of ρ depends somewhat on the placement and shape of the middle voltage contacts. If we have four contacts which are essentially large dots placed along the sample (as is shown in the figure above, and is the case in this lab when we measure VO_2), we will not be able to find ρ exactly. The reason is that our Figure 24 is somewhat idealized in that current may not flow along the sample in such a perfect uniform fashion, being somewhat disrupted by the middle voltage contacts; additionally, we may have a difficulty measuring L precisely, as it is essentially known to the uncertainty in the contact dot size. This is illustrated below in Figure 25.

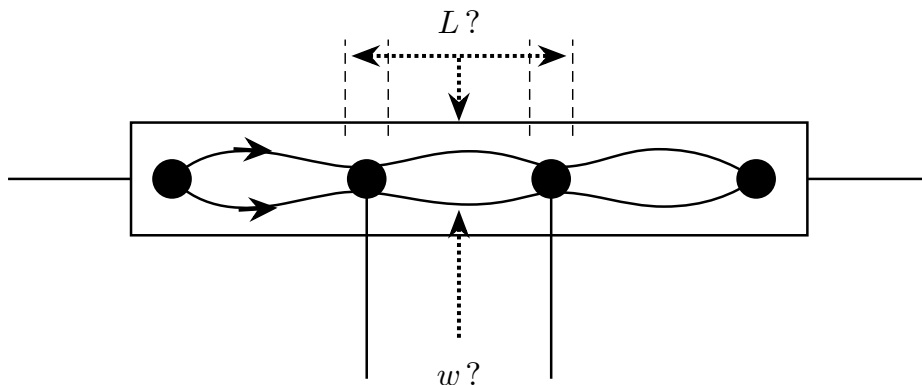


Figure 25: Imperfections associated with placing contacts on the sample surface.

We will ask you to roughly estimate resistivity of the VO_2 sample, despite these imperfections. The thickness δ will be given to you. You may not have perfectly uniform w , but you can estimate average w between voltage leads; and you can roughly measure L between them as well.

We note that one can further improve the four-probe method in order to accurately measure ρ by moving inner voltage leads out of the way of the current and having well defined sample width w . We will not address these improvements here.

B.1 How much voltage-contact resistance can be tolerated?

Finally, how much of a voltage contact resistance can be tolerated before we run into a problem with this measurement? The answer is that resistances in the voltage contacts will become important if they are comparable or exceeding the *impedance of our voltmeter* (in our case, the lock-in amplifier). The input impedance of a voltmeter can be anywhere from $10\text{ M}\Omega$ to a gigaohm; one should consult with the instruction manual of a particular instrument one is using.

To understand why this is so, recall that a real voltmeter *does* need a tiny current for its operation. The last statement applies to both analog and digital voltmeters: whether the current is needed to move the *d'Arsonval galvanometer* in an old analog meter, or to charge a capacitor in a digital one, it is fundamentally needed nonetheless. (The only voltage-measuring instrument which does not need a current is a potentiometer, which is not practical in the type of measurements discussed here.) That is why all voltmeters have finite — not infinite — impedance R_{imp} (internal resistance): they must admit some current. So, a voltmeter is essentially a current-sensing instrument converted to measure a voltage. It fundamentally comprises a current-sensing element (such as a galvanometer G) in series with a very large resistor R_{imp} , as depicted in Figure 26.

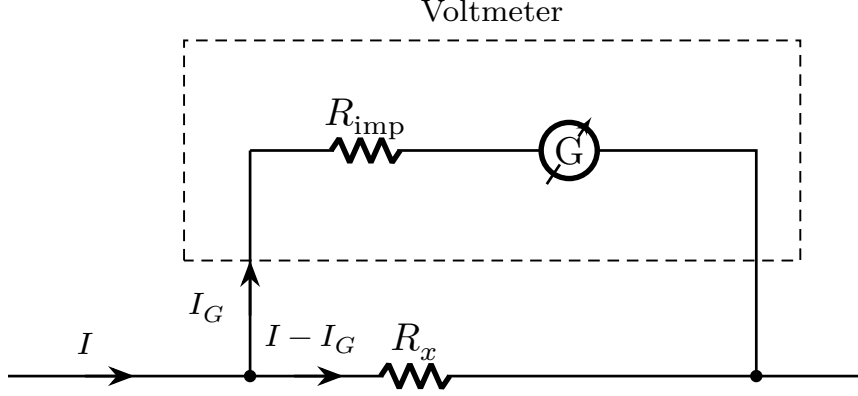


Figure 26: Recalling how a voltmeter works.

We can see that $(I - I_G)R_x = I_G R_{\text{imp}}$, and thus

$$I_G = \frac{(I - I_G)R_x}{R_{\text{imp}}} \approx \frac{IR_x}{R_{\text{imp}}} = \frac{V_x}{R_{\text{imp}}}, \quad (16)$$

where we neglected $I_G \ll I$. The instrument, while measuring I_G , is calibrated to read the voltage across the sample V_x , which is proportional to I_G .

Now suppose that the voltage contact resistances are comparable to R_{imp} . The corresponding electrical circuit showing these unknown contact resistances R_{cont} is shown in Figure 27.

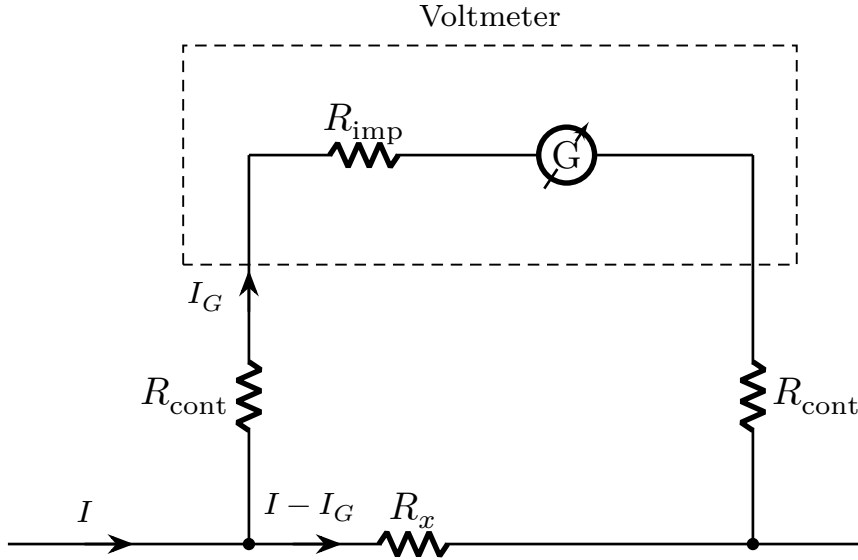


Figure 27: A voltmeter measuring a sample with large contact resistances.

Here we have $(I - I_G)R_x = I_G(R_{\text{imp}} + 2R_{\text{cont}})$, or, again neglecting $I_G \ll I$, $IR_x \approx I_G(R_{\text{imp}} + 2R_{\text{cont}})$, or

$$I_G = \frac{V_x}{R_{\text{imp}} + 2R_{\text{cont}}}. \quad (17)$$

Since R_{cont} is not known, the voltmeter will measure incorrect I_G and thus it will display an incorrect V_x .

Note the difference between this circuit and the one depicted in Figure 23: here a *very small* current I_G flows through contact resistances; there, in a two-terminal measurement, the *full current* I was flowing through these resistances. This is the main advantage of a four-terminal over two-terminal measurement.

Generally it is rather difficult to make good contacts to high-resistivity semiconductors. Contacts to semiconductors and avoidance of contact problems constitute an active field of study, with many publications devoted just to this subject. One characteristic of a good contact is that it is *ohmic*, that is, it has a linear I – V curve, as in Figure 28.

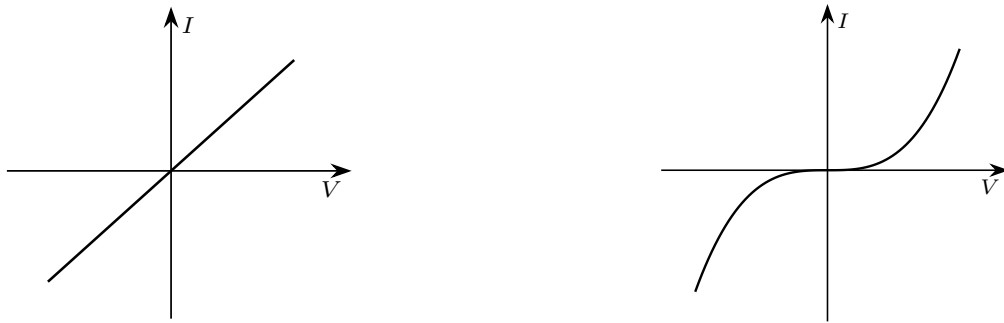


Figure 28: Left: the linear current–voltage curve of a good, ohmic contact. Right: the non-linear curve typical of a high-resistance, “bad” contact.

One often finds that high-resistance “bad” contacts are usually also *non-ohmic*, having a non-linear I – V curve, for example like the one on the right of Figure 28. As one can see, here not only is R_{cont} unknown, but it also is not well-defined, depending on the current I_G that is being passed through it. If in doubt about a contact, measure its I – V characteristic, and stay away from the ones with a non-linear I – V .

Practical advice. When preparing to measure the VO_2 sample, first check its contacts with a simple two-probe ohmmeter. If any of the pairs from the four contacts turn out to have high resistance (say, of the order of several Mohms), report this to the TA or to the Instructor, for you will probably not be able to measure such a sample with the four-probe method.

C Data taking with PyFERRO

PyFERRO replaces the LabVIEW routine `Ferro.vi`. It reads the lock-in, the temperature controller and, optionally, the multimeter, plots the data live and saves them. It only reads the temperature controller: setpoints and ramps are set on the controller itself.

1. Starting PyFERRO. The PyFERRO folder on the desktop contains information about the app. To start PyFERRO, open a terminal in `C:\PyFERRO\pyferro` and run one of the following commands:

- `pixi run start` updates the code to the latest version, then starts PyFERRO with the instruments.
- `pixi run simulate` starts PyFERRO without hardware.
- `pixi run gpib` checks the GPIB connection of the EG&G 5302 (Appendix D); for the SR830, use the Test button below.

The first run downloads Python and all dependencies into the `.pixi\` directory and requires about 1.5 GB of free disk space. Later runs reuse these files. Once setup is complete, PyFERRO can start without an internet connection. If the update step cannot reach the internet, it reports this and starts using the code already present.

2. Instruments tab: Check the lock-in Model (SRS SR830, at GPIB address 8, on this rig), pick the DTECH COM port (FTDI ports listed first), and press **Test lock-in** and **Test controller**. For a hardware run, make sure Simulation mode is off. A green check shows the sensitivity, time constant and reference frequency, and PV/SV, output level, control mode and firmware. A table underneath compares every lock-in setting with the lab manual's (Table 1) and says which key changes any that differ; starting a recording repeats the check and logs any difference as a warning.

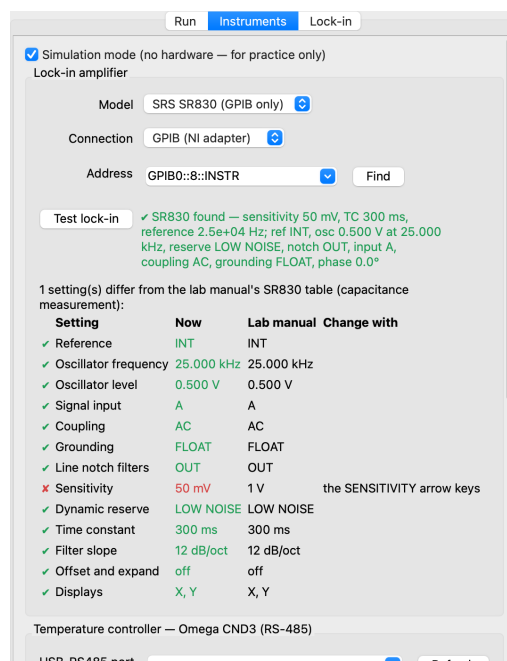


Figure 29: The PyFERRO Instruments tab, shown in simulation mode.

3. Lock-in tab: **Read from lock-in** shows its settings; change them there and press **Apply to lock-in** (it lists the commands and asks first), or press **Lab-manual values** to fill in

the lab manual's table, then Apply. Only what differs is sent, and the settings are read back afterwards. This works during a run too: the change is sent between samples and written into the data file. The phase tuning, and turning an offset or expand on, stay on the front panel.

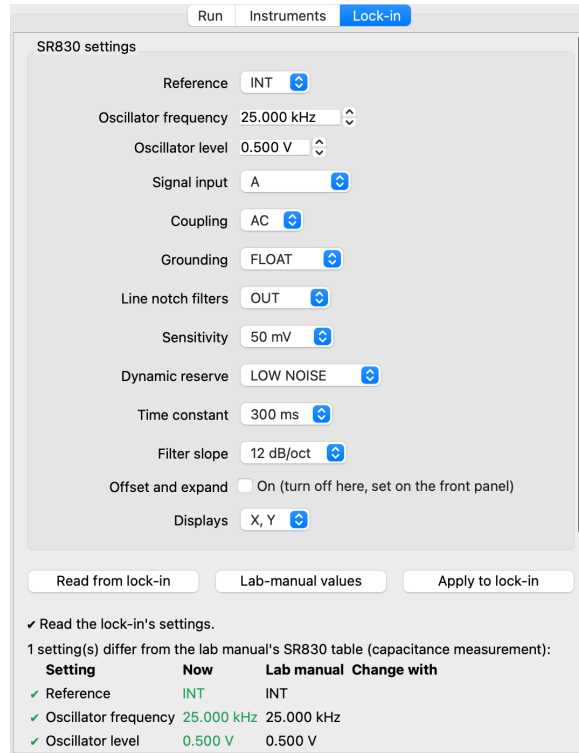


Figure 30: The PyFERRO Lock-in tab for the SRS SR830.

4. Run tab: Enter a sample/run name (becomes the file name), operator, drive details, notes, and the save folder.
5. Monitoring and recording: **Start monitoring** (F5) shows live data without saving. Before starting a new run, press **Clear plots** if you want to clear the displayed plots. Press **Record** (Ctrl+R) to open a file, starting monitoring if needed; pressing it again closes that file, and the next recording opens a new one. **Stop** disconnects.
6. Reading interval and saving: Set the reading interval to at least about five times the lock-in time constant. **Save only if $\Delta T \geq$** reproduces the old LabVIEW behaviour; leave it off to record every reading.
7. Temperature source and multimeter: **Temperature from** (Run tab) selects what goes in column 1 of the file: the CN08D3 controller probe (normal) or the multimeter Pt100. Both are recorded whenever available—PV_C from the controller, T_dmm_C from the multimeter. To use the multimeter, tick **Also read the multimeter** on the Instruments tab and check its model, its GPIB address (GPIB0::24::INSTR for the 34401A, GPIB0::6::INSTR for the Keithley 199 on this rig), and whether the reading is ohms (a Pt100, converted here) or already °C (34401A only).

The tiles above the plots show the sample temperature, setpoint, ramp direction and rate, X , Y and R/θ ; the plots show temperature against time and X and Y against temperature, heating and cooling in different colours. The PyFERRO data file contains (1) the time dependence of the temperature and (2) the temperature dependence of the lock-in output signals X and Y . Copy the data file to a thumb drive for analysis.

Run Instruments Lock-in

Sample / run name simulation

Operator Group / names

Drive / circuit e.g. 2.00 V, 25 kHz, R0 = 600 Ω

Notes Anything needed to reconstruct the measurement (saved in the file header)

Save folder \FerroData Browse... Open

Reading interval 1.0 s

Save only if $\Delta T \geq$ off — save every reading

Temperature from Multimeter Pt100 (GPIB)

Chamber limit warning 160.00 $^{\circ}\text{C}$

☐ Old LabVIEW file format (T, X, Y only; settings in a .json file)

Files are never overwritten: each recording creates *name_YYYYMMDD_HHMMSS.txt* with a header describing the instrument settings.

Figure 31: The PyFERRO Run tab and its data-file settings.

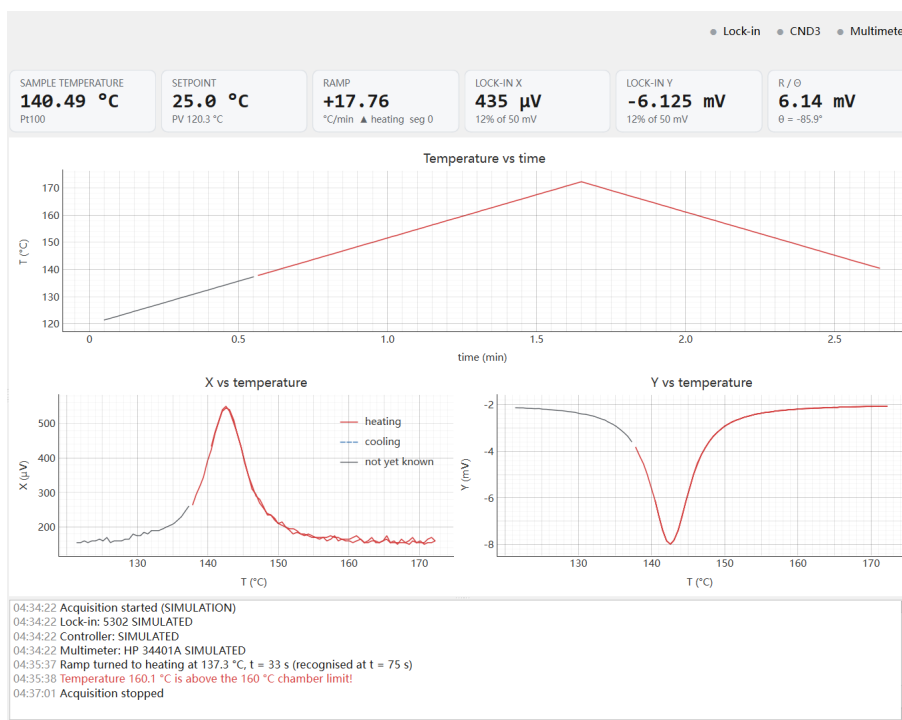


Figure 32: The PyFERRO main window during a simulated run.

A green status light indicates that the instrument is responding; red means that it is not responding (hover over the light for the reason), and grey means that the instrument is unused. A failed reading never stops a run—the value becomes `nan`, a flag is set, and the instrument is reopened after three consecutive failures. The first three temperatures after the controller or multimeter connects are discarded (`nan`, flag 16) so a start-up glitch cannot squash the plot. The temperature tile turns red above the chamber limit (default 160 °C). Settings persist in `%USERPROFILE%\ferro\settings.json`; delete it to reset, or set `FERRO_CONFIG_DIR` for a per-group file.

The data file starts with `#` lines that record the sample name, the lock-in's settings (sensitivity, time constant, reference, oscillator, input, filters and reserve) and the controller's, so the file itself says how the measurement was set up. Anything that changes during the run, such as a sensitivity change or an overload, is written into the file as a `#` line at the moment it happens. The columns are tab-separated; the first three are the temperature (°C) and the lock-in X and Y (V), as in the LabVIEW files, followed by the time, R and θ , the controller's setpoint and temperature, the multimeter temperature, the lock-in full scale, the heating or cooling direction, a segment number and error flags. Missing readings are written as `nan`. The files load directly with `numpy.loadtxt` or with `pandas (comment="#")`, and with most spreadsheet programs.

Write the file name in your lab notes all the same, with anything else you need to reconstruct the circumstances of the measurement.

D The EG&G 5302 lock-in

Before the SR830, the lab used an EG&G Princeton Applied Research model 5302 lock-in, and PyFERRO still supports it (choose it as the lock-in Model on the Instruments tab). If you use the 5302, follow this appendix in place of Section 6.3, and use its own oscillator resistance, $R_0 = 450\ \Omega$, wherever this manual writes R_0 . Its standard settings for the capacitance measurement are:

Setting	Value
Input	preamplifier A input (select PREAMP as signal source)
Sensitivity	1 V initially
Filter	flat
Oscillator level	0.500 V initially
Oscillator frequency	25 kHz
Reference phase	tuned (see below)
Reference	internal
Oscillator out	OSC to preamplifier A (BNC cable)
Time constant	500 ms
Dynamic reserve	HI STAB
Output display	X-YV

Setting the EG&G model 5302. Table 2 gives the same settings in more detail: where to find each one on the front panel, and the equivalent GPIB command. The same command sent without a number makes the lock-in reply with the current value, which is how PyFERRO reads these settings back and records them in the data-file header.

Oscillator resistance and level. The 5302’s oscillator output resistance is $R_0 = 450\ \Omega$ (5302 manual, section 4.6), not the $50\ \Omega$ of the SR830: use $450\ \Omega$ in $C = \delta/(\omega R_0)$. Its least sensitive range is 1 V, and the in-phase signal is close to the oscillator level V_0 , so V_0 must stay below 1 V; at $V_0 = 0.500\ \text{V}$ it is half of full scale. Nothing is lost, since C depends only on the ratio $\delta = V_{\text{meas}}/V_0$.

The 5302 is a dual-phase lock-in: it shows the in-phase signal as X and the out-of-phase signal as Y at the same time. Tune the phase with the sample disconnected until Y is as close to zero as possible (Auto-Phase does this), then connect the sample and read V_{meas} from Y .

The practice run with a piece of BNC cable (Section 6.3) goes as follows with the 5302.

First, check the signal generated by the lock-in amplifier’s internal oscillator. Connect one end of a BNC cable to the “OSC” output and the other end to input A on the preamplifier. Select “PREAMP” as the signal source. On the main display, select “PHASE” and turn the setting knob to set the reference phase to 0.0° . Check that the analog meter and the in-phase component X show a clear signal. If most of the signal appears in X , its reading should be close to the oscillator setting of $0.500\ \text{V}$.

We are interested in the out-of-phase component Y . Ideally, Y should be zero after the phase has been correctly adjusted. However, the electronics inside the lock-in amplifier can introduce a small phase error, so Y may initially have a small but nonzero value.

The lock-in amplifier provides an automatic way to correct this phase error. Press the red “AUTO” button, followed by the left “PHASE” button, to run Auto-Phase. Wait a few seconds for the readings to settle. The in-phase component X should then be maximized, while the out-of-phase component Y should be close to zero.

Upon closer inspection, you may find that Auto-Phase has left a small out-of-phase signal. If Y is too small to read clearly, keep the sensitivity at 1 V and press the “+90” phase button once to bring this small signal into the X channel. Then press “FUNCT” followed by “SEN down/EXPAND” to magnify X by a factor of ten. Watch X and adjust the reference phase manually until

Setting	Value	How to set it on the front panel	GPIB
Signal input	preamplifier, input A	Cable to input A of the plug-in preamplifier. The display shows PRE at the upper right; if it shows DIR, select PREAMP as the signal source on the SIGNAL SETUP screen.	PREAMP 1
Reference	INT	Press REF, the rightmost of the five keys under the display, until INT shows in the upper-left window. The INT light above the meter comes on.	IE 0
Oscillator frequency	25.000 kHz	Press the keys to the right of the large setting knob until OSC F is in the middle of the lower-right window, then turn the knob. The KNOB ACTIVE light must be on; if it is off, press either key once.	OF 2500 7
Oscillator level	0.500 V	As for the frequency, with OSC V selected.	OA 5000 1
Sensitivity	1 V	The SEN keys, at the top of the two key columns to the right of the display: the left one moves to a larger full scale (less sensitive).	SEN 21
Expand	off	When it is on, EX shows at the upper right. Press FUNCT and then the right SEN key (SEN/EXPAND) to toggle it.	EX 0
Time constant	500 ms (or 200 ms)	The TC keys, below the SEN keys. The 5302 has no 300 ms.	XTC 8 (XTC 7)
Filter	FLAT	Press FILT, the second key from the right under the display, until FLAT shows at the upper left.	FLT 0
Dynamic reserve	HI STAB	The DYNRES/LOCAL key, to the left of the display near its bottom edge.	DR 1
Reference phase	$Y = 0$ with the sample disconnected	Press AUTO and then the left PHASE key (Auto-Phase); repeat once the outputs have settled. For fine tuning, select PHASE with the keys by the setting knob and turn the knob.	AQN
Display	X-YV	The DISP keys, to the right of the display.	—

Table 2: EG&G model 5302 settings for the capacitance measurement. The GPIB commands are from chapter 9 of the 5302 manual; OF 2500 7 means 2500 in the 10 to 100 kHz range and OA 5000 1 means 5000 steps of 0.1 mV in the 25 to 500 mV range.

it is as close to zero as possible. Press “FUNCT” followed by “SEN down/EXPAND” again to turn off EXPAND, then select “PHASE” and subtract 90° from its current setting to restore the original X-Y orientation. If the red “OUTPUT OVLD” indicator lights while EXPAND is on, turn EXPAND off before continuing. Do not reduce the 1 V sensitivity setting while the oscillator is connected directly to the input.

Now that the out-of-phase signal has been brought as close to zero as possible, keep the sensitivity at 1 V. Insert a BNC “T” connector into the line and connect another BNC cable. The added cable will produce a nonzero out-of-phase signal V_{meas} , which you can read from Y. Calculate its capacitance using $C = |V_{\text{meas}}|/(V_0\omega R_0)$, where $V_0 = 0.500$ V, $\omega = 1.57 \times 10^5$ s $^{-1}$, and $R_0 = 450 \Omega$. Measure the length of the cable, then repeat the measurement with a longer cable.

E Physics of phase transitions

E.1 Introduction to phase transitions

An interesting, important and useful class of phenomena goes under the name *phase transition* [5, 6]. Phase transitions are widespread in nature, and a number of them find useful applications



Figure 33: The HP 34401A multimeter (top) and the EG&G Princeton Applied Research model 5302 lock-in amplifier (bottom).

in technology. In this lab we will primarily concentrate on phase transitions which manifest themselves in electrical properties of solids.

Microscopic theories of different phase transitions involve diverse phenomena specific to these transitions. Yet there are some similarities which allow their classification and a (somewhat) unified description. The phenomenological theory — i.e. a non-microscopic theory, one which mathematically describes given phenomena without asking detailed questions about their mechanism — of phase transitions was developed largely by a famous Soviet physicist Lev Landau around 1937 [5, 6, 7, 8]. It is based on thermodynamic ideas and on symmetry. While the detailed description of this theory exceeds the level at which we can present the subject here, some of the Landau ideas are actually quite simple and transparent, and we will present them here, specifically in Section E.3.2.

We will start with a general principle of phase transition classification. Then, the following sections will detail two phase transitions: the metal–insulator transition [9, 10, 11] and the ferroelectric transition [12, 13, 14] that you will study in this lab.

E.2 Classification of phase transitions

A phase transition is usually defined as an abrupt change in physical properties upon a *continuous change in external parameters*. Examples of external parameters are temperature, pressure, magnetic and electric fields. Of these temperature is the most common cause of a phase transition, and in this lab you will be studying such temperature-induced transitions. The value of the temperature, pressure, or other physical quantity at which a phase transition occurs is called a *transition point*. If it is temperature, it is usually denoted as T_C .

A distinction is made between two general classes, or, as they are commonly called, *orders* of phase transitions. This distinction follows pioneering classification of phase transitions proposed in the 1920s by Paul Ehrenfest. Why do we speak of a transition order? The reason is that Ehrenfest

classified phase transitions based on the behavior of the thermodynamic free energy [5, 6] (see below) as a function of other thermodynamic variables. Under this scheme, phase transitions were labeled by the lowest derivative of the free energy that is discontinuous at the transition. Thus, according to this scheme, *first-order* phase transitions are defined as transitions exhibiting a discontinuity in the first derivative of the free energy with respect to some thermodynamic variable, and *second-order* phase transitions are continuous in the first derivative of the free energy with respect to some thermodynamic variable, but exhibit discontinuity in a *second derivative*. Free energy itself is always continuous at the transition point, which is required by conservation of energy.

A modern classification is only slightly different from Ehrenfest's. The meaning of transition's order in terms of discontinuous free energy derivatives remains, while the distinction between first and second-order transitions is formulated as follows: in a first-order phase transition, specified amount of heat, called the *heat of transformation*, or *latent heat*, is released or absorbed per unit mass. In a second-order phase transition, there is no latent heat — heat is neither released nor absorbed — but instead some physical quantity, called an *order parameter* appears; order parameter is equal to zero on one side of a transition point (in other words, it does not exist in one of the phases), and it gradually increases from zero on the other side. It was Lev Landau, in his above-mentioned theory of second-order phase transitions [5, 6, 7, 8], who first raised the order parameter to prominence. We will say more on this subject later.

Examples of first-order phase transitions include many of the most familiar ones, such as the water–ice transition, or the water–vapor transition; more generally — all the solid/liquid/gas transitions are of first order. An interesting and practically important first-order phase transition manifests itself in the electrical and optical properties of a solid: it is a *metal–insulator transition* [9, 10, 11], or, as it is sometimes called, a semiconductor-to-metal transition. It can be found in a fairly large number of solid compounds (but never in pure elements), including vanadium dioxide (VO_2), the material you will be studying in this lab.

Second-order phase transitions also include many of the transitions of great interest in physics and its applications. A familiar and important example is a *ferromagnetic transition* [5, 6, 13, 15] at a Curie point: above the Curie temperature (in iron it is $T_C = 770^\circ\text{C}$) the material is in a paramagnetic state, with randomly oriented magnetic moments; at the Curie point the material forms domains with magnetic moments all aligned in the same direction within a domain. This phenomenon is called *spontaneous magnetization*. The order parameter we mentioned above in this case is the net magnetization inside a domain; it is zero at T_C and increases from zero below T_C . We note also that zero value of order parameter corresponds to a completely disordered state (random orientation of magnetic moments), while its finite value corresponds to an appearance of order (common direction for all the magnetic moments in a domain). At the same time, the material in paramagnetic state is more symmetrical than in the state of spontaneous magnetization: rotational symmetry of the macroscopic collection of randomly-oriented microscopic magnetic moments in a paramagnetic state is broken when material transitions into the oriented state below T_C . Thus such a phase transition is intimately related to symmetry, the fact which was fully explored by Landau in his theory.

We note here that the ferroelectric transition in barium titanate (BaTiO_3) which you will study in this lab is similar to the ferromagnetic one, with electric dipole moments and electrical polarization taking the place of their magnetic counterparts (we will discuss the differences in Section E.7). We also note that ideas of spontaneous symmetry breaking migrated from solid state physics into cosmology and elementary particle theory, and are presently of singular importance in these fields, illuminating a large number of fundamental puzzles; for example, the celebrated Higgs mechanism responsible for the appearance of mass was first proposed by a famous solid state theorist P. W. Anderson in 1962, who at the time worked on the aspects of theory of superconductivity.

Another example of a second-order phase transition is a *superconducting transition* [15], which

takes place below a transition temperature (also called the critical temperature) T_C : at $T > T_C$ one observes a normal, resistive metal; at $T = T_C$ the electrical resistance abruptly (in pure metals) disappears and a normal metal turns into a superconductor with zero resistance; at lower temperatures $T < T_C$ the resistance is identically zero. The order parameter in this case is the *superconducting gap*: an energy gap which is equal to zero at and above T_C and opens up at lower temperatures. It is the gap in the electronic spectrum, of the order of a meV, separating paired electronic state (so-called Cooper pairs) from the normal electron spectrum. Note that superconducting gap is quite different from the gap in a semiconductor or insulator electronic spectrum, of which we will speak in Section E.4. The entropy of a superconducting state is lower than that of a normal metal state, indicating that electron pairing provides for a more ordered state.

E.3 Thermodynamic approach to phase transitions

In discussing phase transitions, it is instructive to take a thermodynamic approach, considering Helmholtz free energy

$$F = U - TS \quad (18)$$

of a material at the two sides of a phase transition. Here U is internal energy, T is temperature, and S is entropy. Sometimes instead of Helmholtz free energy F , Gibbs free energy G is being considered; it is equal to $F + pV$, where p is pressure and V is volume. In a phase transition there are generally two phases, one being stable on one side of the transition, another on the other side. Thermodynamically, a phase transition occurs at the external parameter value (temperature, or other) at which *the free energies of the two phases become equal* (which guarantees, as we said above, that the free energy is continuous across the transition). One can make some conclusions about the two general types of phase transitions by making some general observations concerning formula (18).

E.3.1 First-order transitions from the thermodynamic point of view; mixed-phase states

We recall the First Law of Thermodynamics which states that the heat Q absorbed by a system is equal to the change in internal energy of a system ΔU plus the work done by a system W : $Q = \Delta U + W$. If the system does not produce any work (i.e. if it does not expand or contract against some outside pressure), then $Q = \Delta U$. This is generally the case in solid-state transitions of interest to us here. This tells us that first-order phase transition, which, as we said above, is accompanied by the appearance of latent heat Q , is possible only if internal energy of the system changes. This energy change must be electronic or ionic in nature: either the electron energy spectrum, or the electrostatic energy of the ionic lattice, or both, must undergo a change at the transition point. For example, there may be an energy gap E_g opening up on one side of a phase transition, as is the case in VO_2 . Such a rearrangement in the energy spectrum of a solid will be abrupt: in our example, the semiconductor energy gap in VO_2 exists on one side of the transition and does not exist on the other; in other words, it does not gradually open up starting from zero at the transition, as the superconducting gap. The latter, as we said, serves as an order parameter in that *second-order* transition. Semiconducting gap, in contrast, could not be considered an order parameter.

Electronic energy spectra in VO_2 on the two sides of a transition are shown schematically in Figure 34.

Further, we see that in formula (18) there are two terms: U and TS . We already decided that in a first-order transition U must change. But we also know that the free energies of the two phases must be equal at T_C . Therefore the second term, TS , must also change at the transition in order to compensate for ΔU . Indeed, S , the entropy, must change if there is latent heat Q , being related to it by $\Delta S = Q/T$. Change of entropy signifies appearance or disappearance

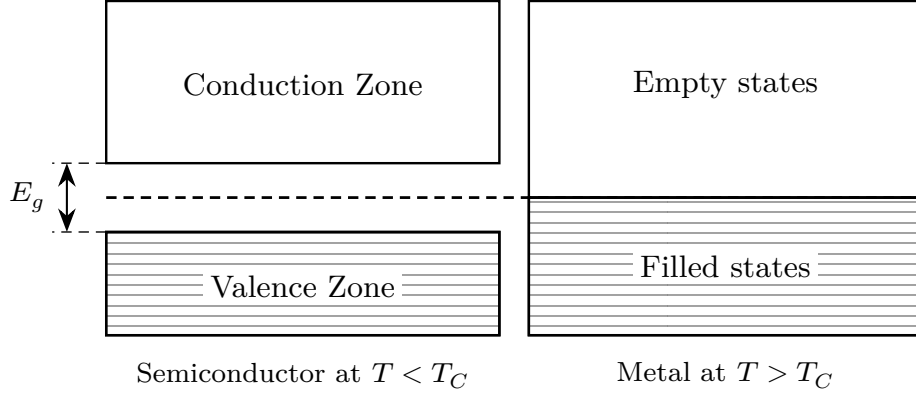


Figure 34: A diagram schematically representing the electron energy spectrum in a first-order phase transition in which an energy gap opens below T_C ; such a transition is observed in VO_2 .

of order (here the word *order* is used in a context of order-disorder and not in a context of derivative's order). For example, in VO_2 one finds different lattice symmetry in the two phases: one phase is tetragonal, another monoclinic (these terms refer to types of crystalline atomic arrangements). And indeed one finds a change in entropy which comes predominantly from this lattice rearrangement (see below a section on VO_2). The lattice constant of a unit cell and material's density are also changing in the first-order transition. Note that this is different for the second-order transition, where *order* does not change at T_C , but changes gradually away from T_C .

Let us continue with general considerations concerning a first-order transition. We are considering two phases, phase 1 and phase 2, on the two sides of a transition. Suppose we are crossing the transition temperature T_C going from phase 1 to phase 2. The free energy, as we said above, is a continuous function across the transition, so that $F_1(T_C) = F_2(T_C)$. However, $F(T)$ is not constant; each phase has lower F on its “rightful” side of a transition. This is schematically shown in Figure 35, where phase 1 is stable at $T < T_C$ and phase 2 is stable at $T > T_C$.

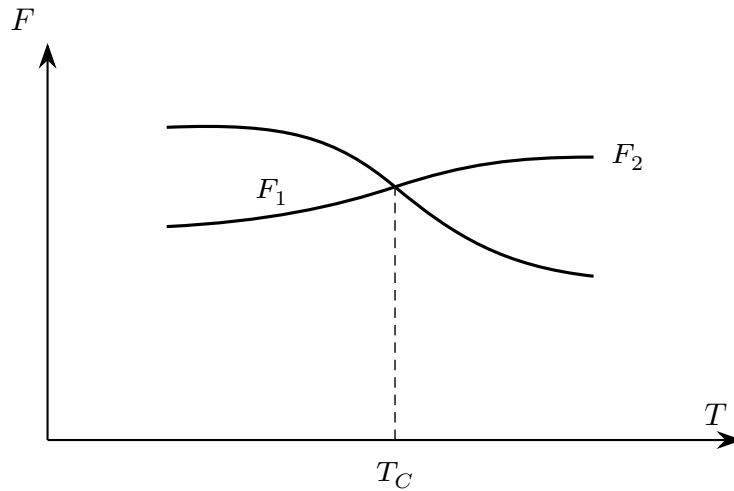


Figure 35: Schematic behavior of the free energy of the two phases across the phase transition.

If a given phase — for example, phase 1, is energetically favored (having lower F), one might expect it to appear right away at crossing the T_C . But consider in more detail the mechanism

of its appearance: what has to happen for it to appear? Some atomic rearrangement and/or lattice symmetry change will be required; electronic and ionic energy spectra will have to shift; latent heat will have to be released or absorbed. Very importantly, unless the whole volume of the material uniformly transitions all at once, domains of phase 1 may form, surrounded by the phase 2, and, as a result, there appears the *boundary* between the two phases. This boundary may have energy, either positive or negative. In case boundary energy is positive, the boundary is said to “cost” energy. Appearance of a boundary thus is changing somewhat the energy balance dictated by the lower free energy of phase 1. There may appear a *local energy minimum*, which depends on the local configuration of the phase domains; this local minimum may preclude the system from falling into the *global minimum* dictated by the lower F_1 value. Thus the phase transition may not happen exactly at T_C , and it may develop in a complex, non-uniform fashion. During this process, the temperature of the system will stay constant as heat is added: the system is in a *mixed-phase regime* in which some parts of the system have completed the transition and others have not.

Such phenomena are well known in first-order phase transitions between solids and liquids and liquids and vapor: a liquid may be heated to a temperature above the boiling point or *supercooled* to below the freezing point. Something similar happens in a first-order phase transition in a solid. First-order transitions are often characterized by the formation of *metastable* states. Instead of taking place exactly at T_C , the transition going up in temperature takes place at $T_C + T^*$, and going down in temperature at $T_C - T^*$. Here T^* is called the *coercive temperature*. This leads to the phenomenon of *hysteresis*. If we fix a temperature T somewhere in the interval $T_C - T^* < T < T_C + T^*$, a metastable state, often in the form of domains, may survive for long (sometimes indefinitely long) periods of time on the “wrong” side of the transition. In Section E.4.5 you will find a qualitative discussion of the reason for the appearance of a coercive temperature and hysteresis.

In this lab you will study one such first-order phase transition, namely the semiconductor-to-metal transition in vanadium dioxide (VO_2), which also exhibits a mixed-phase regime. Additionally, the semiconductor phase of VO_2 has high resistivity, while its metallic phase has low resistivity. Thus, as we explore the mixed-state with electrical resistance measurement, we are encountering another interesting type of a transition here, called *percolation* [16] transition. We will briefly discuss these phenomena when talking about VO_2 in Section E.4. In this lab, when performing measurements on VO_2 , you will be able to see how, with increasing temperature, resistivity gradually changes from high to low, dramatically changing by several orders of magnitude in the process. The material goes through the metal-insulator transition, a two-phase mixture occurs, and, as a result, it also goes through the percolation transition.

E.3.2 Second-order transitions from the thermodynamic point of view; elements of Landau theory

In a second-order transition there is no latent heat, that is, $Q = 0$ and thus, from the First Law of Thermodynamics, $\Delta U = 0$. Thus an energy spectrum is not expected to modify abruptly at the transition point in a second-order transition. At the same time, the general condition $F_1 = F_2$ at $T = T_C$ must be satisfied. Looking at formula (18) for F , this tells us immediately that entropy (order) cannot change at T_C , but has to change gradually as the temperature deviates from T_C . In other words, the order parameter must be zero at T_C , and will increase away from T_C , as we said above. This, among other considerations, allowed Landau to construct his theory of the second-order phase transitions. Namely, the fact that order parameter is small in the vicinity of a transition point allowed him to expand free energy F as a Taylor series in order parameter, which we will call m , and then, making general physical arguments, to extract significant physical results from this expansion.

Let us follow his lines of thought, at least in part, here. The Taylor expansion of a continuous function $F(T, m)$ in powers of a small parameter m around point $m = 0$, written to the 4-th

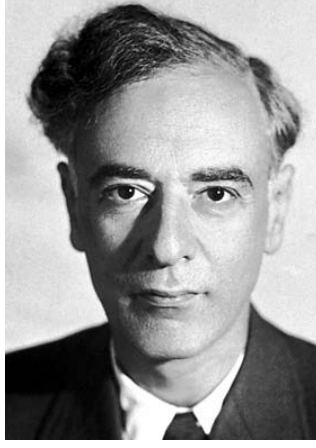


Figure 36: L. D. Landau (1908–1968), Nobel Prize of 1962.

order of m , is

$$F(T, m) = F_0(T, 0) + \frac{F'(T, 0)}{1!}m + \frac{F''(T, 0)}{2!}m^2 + \frac{F'''(T, 0)}{3!}m^3 + \frac{F''''(T, 0)}{4!}m^4 + \dots \quad (19)$$

Here primes indicate derivatives with respect to m ; for example $F'' = d^2F/dm^2$. So far this is just a mathematical formula based on continuity of F and smallness of m . Landau then argued that F should be invariant with respect to a change in sign in m : changing m to $-m$ should keep the energy unchanged. Indeed, in the absence of external field, if m is magnetization M or polarization P , it seems obvious. This is a symmetry argument we mentioned before. This implies that our Taylor expansion must not contain odd powers of m . Omitting odd terms with m and m^3 , and calling the even-term coefficients $F''(T, 0)/2! = \alpha(T)$ and $F''''(T, 0)/4! = \frac{1}{2}\beta(T)$ (the factor $1/2$ introduced here following a convention found in most of the literature), we obtain

$$F(T, m) = F_0(T) + \alpha(T)m^2 + \frac{1}{2}\beta(T)m^4 + \dots \quad (20)$$

We want to relate behavior of $F(T, m)$ to the observed physical phenomena in the vicinity of a phase transition. *A minimum of F with respect to an order parameter m physically implies stable (observable) configuration of our system:* the value of magnetization or polarization observed in a real system will correspond to the value of m at which F has a minimum. We will find that these minima are different for $T > T_C$ and for $T < T_C$, on the two sides of a phase transition. After finding the minima with respect to m , we may be able to see how F depends on m in the vicinity of these minima. Once we decide on the form of $F(T, m)$ for fixed T , we will be in the position to discuss T -dependences. Further, by taking derivatives of $F(T, m)$, one can construct other quantities, such as entropy, specific heat, etc. Such is our “program” of how to explore physical consequences of expansion (20).

To find minima of F with respect to the order parameter m , we require that $dF/dm = 0$ and $d^2F/dm^2 > 0$: the usual conditions we always use when finding a minimum of a function. The first one, $dF/dm = 0$, guarantees that we have zero slope, i.e. that we are at the extrema of a function; the second condition, $d^2F/dm^2 > 0$, guarantees that this extremum is actually a minimum and not a maximum or the inflection point.

Differentiating (20) we find

$$\frac{dF}{dm} = 2\alpha(T)m + 2\beta(T)m^3 = 2m [\alpha(T) + \beta(T)m^2] = 0, \quad (21)$$

which gives us three possible solutions for m :

$$m_1 = 0 \quad \text{and} \quad m_{2,3} = \pm \left[-\frac{\alpha(T)}{\beta(T)} \right]^{1/2}. \quad (22)$$

The second condition for a minimum is $d^2F/dm^2 = 2\alpha(T) + 6\beta(T)m^2 > 0$. The simplest way to satisfy this last inequality is to require that at $T > T_C$, $\alpha(T) = a(T - T_C)$ with coefficient $a > 0$, and that $\beta(T)$ is a positive constant $b > 0$. We note that if we adopt this form for $\alpha(T)$, then $\alpha(T_C) = 0$.

So, now the expansion looks like this:

$$F = F_0(T) + a(T - T_C)m^2 + \frac{1}{2}bm^4 + \dots \quad (23)$$

And F has minima at

$$m_1 = 0 \quad \text{and} \quad m_{2,3} = \pm \left[\frac{a}{b} (T_C - T) \right]^{1/2}. \quad (24)$$

Clearly the expression for $m_{2,3}$, by virtue of containing a square root, works only for $T < T_C$ (temperatures below the transition); thus $m_1 = 0$ corresponds to $T > T_C$ (temperatures above the transition). As can be seen from (24), at the transition temperature $m = 0$ in both cases.

Now we will try to explore the physical consequences. Recalling the physical meaning of m , we see that equilibrium magnetization (or polarization) is zero at $T > T_C$: this is clearly the description of a paramagnetic (paraelectric) state, in which magnetization from randomly-oriented magnetic moments equals to zero. According to the expansion, at $m = 0$, $F = F_0(T)$. At $m \neq 0$, and with positive $T - T_C$ (above the transition), expansion (23) predicts that F will rise from the lowest energy level F_0 as $F - F_0 = Am^2 + Bm^4$, where A and B are positive coefficients. This is shown schematically in Figure 37(a), where we took F_0 to be zero.

The minima corresponding to $m_{2,3}$ describe the behavior of F below the phase transition: our theory predicts that there will be two such minima, symmetrical in m with respect to $m = 0$. Between these two minima, at $m = 0$, in virtue of (23), we have $F = F_0$ (or zero, since we have chosen zero at F_0). This is depicted in Figure 37(b). Physically this means that our system will have two possible stable states with non-zero magnetization or polarization. Thus in a ferromagnetic (or ferroelectric) state magnetization (polarization) will have a *finite* value. This finite value may be positive or negative; it can reverse, still being consistent with our theory.

The values of F at these minima can be calculated by substituting m_1 and $m_{2,3}$ into (23). They will depend on temperature T , which is now less than T_C . We see from (23) that the closer is T to T_C , the smaller will be the negative term $a(T - T_C)m^2$ in the expansion, and thus the shallower will be the two minima below the transition. At the same time, with decreasing $T - T_C$, the two minima will move in towards the $m = 0$ value. In the limit $T \rightarrow T_C$ the two minima will disappear, transforming the double-minima curve into a single minimum curve (a). At $T = T_C$, because of the m^4 term, the curve will have a significant flat region in the vicinity of $m = 0$. Physically this must imply that, near the transition point, magnetization (polarization) in a ferromagnetic (ferroelectric) state can be very easily reversed. In other words, the material at a temperature just below the transition point must be very easily polarizable. This suggests that *dielectric constant* ε will be high. We will see that indeed this interesting conclusion is experimentally verified, and you will be able to see this very clearly in your measurements of BaTiO₃.

So far our conclusions based on the analysis of expansion (23) have been very satisfying: we were able to explain some of the most prominent features of a second-order phase transition. We could continue and explore the physical consequences for the derivatives of F , which will give us entropy and specific heat. We will find [5, 6] that entropy S , while being continuous at the transition, decreases below the transition, where there is more order. We will also find that specific heat experiences a jump at the transition. Such a discontinuity in specific heat at the second-order phase transition is well documented experimentally.

One can add external field (magnetic field H , or electric field E) to this analysis. In doing so one finds [5, 6] that susceptibility χ , which is defined as $\chi_M = dM/dH$ for the magnetic case and $\chi_E = dP/dE$ for the electric case [5, 6, 15], obeys the Curie–Weiss law (here written for the

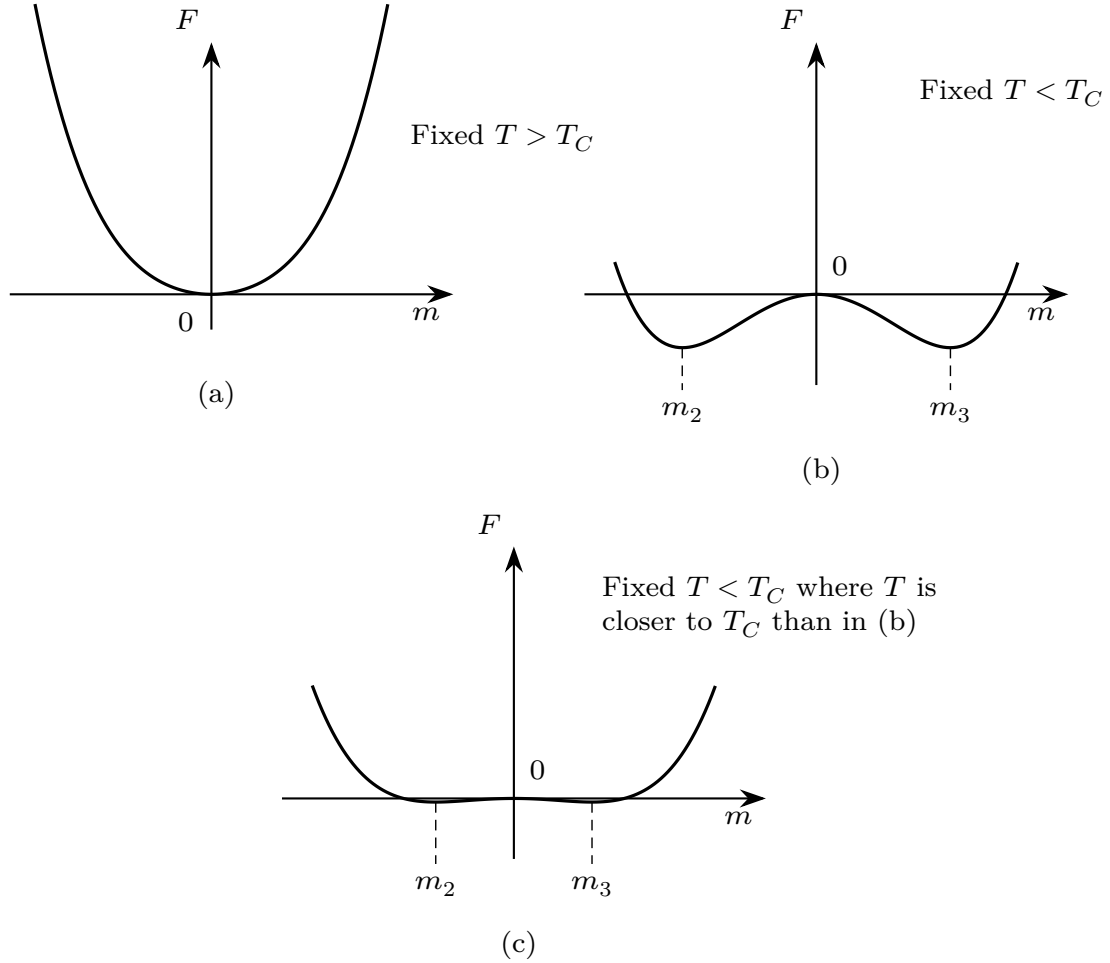


Figure 37: Schematic representation of the dependence of the free energy F on the order parameter m at three fixed temperatures.

ferroelectric case):

$$\chi_E = \frac{1}{2a(T_C - T)} \quad \text{for } T < T_C, \quad \chi_E = \frac{1}{4a(T - T_C)} \quad \text{for } T > T_C. \quad (25)$$

This shows that Landau theory predicts that susceptibility becomes infinite at $T = T_C$. We also see that it approaches $T = T_C$ in an asymmetric way (constants being different by a factor of 2). In an experiment you will perform in this lab you will see that this is indeed so, except that theoretical infinity is substituted by a sharp peak.

We recall Curie–Weiss law [5, 6, 15] which states that $\chi = C_{CW}/(T - T_C)$ for the paramagnetic (paraelectric) phase, i.e. for $T > T_C$. Here C_{CW} is a Curie–Weiss constant (often called Curie constant). We see that Landau theory identifies coefficient $1/(4a)$ with C_{CW} .

E.4 Description of the phase transitions and materials you will study

E.4.1 Metals, semiconductors and insulators, and their conductivity

With respect to an ability to transmit electrical charge, materials can be divided into conductors and insulators; this was recognized from the eighteenth century. Ohm's law, resistance and resistivity were defined and clarified in the nineteenth and the beginning of the twentieth century. Today, with respect to electrical resistivity ρ (or its inverse $1/\rho = \sigma$, electrical conductivity)

we distinguish not two but three classes of materials: metals (good conductors), dielectrics or insulators (very bad conductors), and semiconductors, which are in some sense intermediate between the first two classes. We also distinguish them by different ways in which conductivity depends on temperature (below). True understanding of what is at the core of these different classes came only in the 1930s, when physicists understood electronic spectra (band structure) of these materials [13, 15]. This required creation of quantum mechanics, formulation of Pauli's exclusion principle, and replacement of classical statistics with Fermi–Dirac statistics. In Figure 38 we show a schematic representation of the electronic bands of a metal, a semiconductor and an insulator.

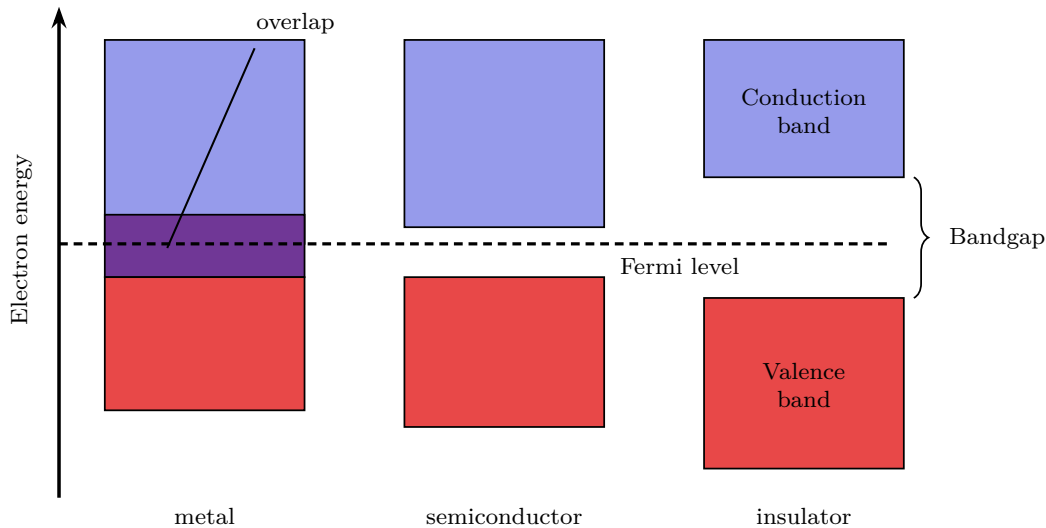


Figure 38: Electronic spectra (band structure) of metals, semiconductors and insulators.

The main distinction between them, which influences electrical conductivity and its temperature dependence, lies in the presence (in semiconductors and insulators) or absence (in metals) of an energy gap in the electronic spectrum. The value of the energy gap is called E_g . The highest energy band which contains states which are completely filled with electrons is called a *valence band*; an empty band, separated from it by an energy gap, is called a *conduction band*. Electronic states of the type that may conduct current (the so-called extended or Bloch states) do not exist in the energy interval inside the gap.

The difference between semiconductors and insulators is quantitative rather than qualitative: they both have a gap in their spectra, but the gap in a semiconductor is smaller than that in an insulator. The rough borderline between these two classes lies around $E_g = 2$ to 3 eV. One can say that semiconductor materials are small band-gap insulators. It should be said that, in addition to the presence of a relatively small gap, another defining property of a semiconductor material is that it can be doped with impurities that alter its electronic properties in a controllable way. Here we will be concerned mainly with the so-called *intrinsic* semiconductors, which do not have impurity levels inside the gap. The most famous semiconductors today are Si and GaAs (and their extended families); but there are hundreds of known semiconductors, and many of them find useful applications. VO_2 you will study in this lab is one of them.

Let us understand why conduction will be so very different depending on the presence or absence of a gap in the electronic energy spectrum. You will need to understand this, as we will ask you to find the value of the gap in VO_2 from your data.

In a conductor, application of an electric field E leads to acceleration of charge carriers (such as electrons) with mass m and charge e : $a = eE/m$. This acceleration acts for some average time between collisions τ , leading to the appearance of a drift velocity $v_d = a\tau = (eE/m)\tau$. The

current density is then

$$j = env_d = \frac{e^2 n \tau}{m} E = \sigma E \quad (\text{Ohm's law}). \quad (26)$$

Here n is the number of carriers per unit volume (carrier density), and σ is electrical conductivity which characterizes the material. Its inverse is resistivity ρ .

It is clear from these familiar considerations (constituting part of what is known as Drude model [13, 15] of electrical conduction) that an electron must be able to *change* its energy in order to participate in electrical current; indeed, it acquires extra kinetic energy $mv_d^2/2$ from the field; this extra energy will eventually be lost in inelastic collisions of an electron with phonons (lattice vibrations), other electrons, etc. We know that eventually this extra energy is deposited into the conducting material, raising its temperature: the familiar $I^2 R$ power dissipated in a conductor comes from the field (the battery) and ends up as heat in the conducting material. Note that energy increase $mv_d^2/2$ is, generally, quite small compared to the total kinetic energy of an electron at the highest filled energy level.

So, to repeat: an electron needs to change its energy in order to accelerate under the influence of an electric field; this is the prerequisite for electrical conduction. But, to change its energy, it has to find un-occupied energy levels at that higher energy. This is possible in a metal, where un-occupied states lie directly above the occupied ones (at zero temperature we call the highest energy of occupied states *Fermi level*; it is shown as a dashed line in Figure 38; there are empty states right above it). In a semiconductor or insulator, in view of the existence of a forbidden energy interval (energy gap) this is impossible: there are no available states above the top edge of the valence band (see Figure 38). On further reflection, however, we see that this last statement is strictly true only at $T = 0$. Indeed, at a finite temperature T electrons may be excited across the gap. Relative probabilities of such processes can be found from Fermi–Dirac distribution functions; approximate (but very good) estimates can be done, however, using classical Boltzmann factor $\exp(-E_g/kT)$. Such an excitation process creates an electron near the bottom of the conduction band above the gap, and a *hole* (which is equivalent to a positively charged electron [6]) near the top of the valence band. These charge carriers (electron and hole) are called *mobile carriers*; they can now participate in conducting the current, because they can find un-occupied energy levels very close to their respective energies.

The number of mobile electrons and holes depends on the values of T and E_g , and on how many electrons were available for excitation. In good insulators, i.e. materials with $E_g > 2$ eV, at practical temperatures (limited, generally, by the melting point; often we are interested in conductivity at $T = 300$ K), their numbers are negligibly small. In other words, the resistivity of an insulator is many orders of magnitude higher than in a semiconductor or in a metal. For example, resistivity of a very good insulator amber (the one that ancient Greeks were rubbing to produce static electricity) at 300 K is $5 \times 10^{14} \Omega \text{ m}$, while resistivity of pure Si is $2300 \Omega \text{ m}$, a difference of eleven orders of magnitude! Metals are still orders of magnitude more conductive than semiconductors: resistivity of copper at 300 K is only $1.7 \times 10^{-8} \Omega \text{ m}$. Thus metals and good insulators differ by staggering 22 orders of magnitude (!) in resistivity (or conductivity).

In addition to very different values of ρ , its temperature dependencies are very different as well: in good metals resistivity increases with T ; in copper and many other metals this increase is linear at sufficiently high T . In contrast, resistivity in a semiconductor depends on temperature as

$$\rho = \rho_0 \exp\left(\frac{E_g}{2kT}\right). \quad (27)$$

Here the factor of 2 comes from the fact that activation energy is counted from the Fermi level, or, as it is called in semiconductors, from chemical potential, which, in intrinsic semiconductors, lies in the middle of the gap [13, 15]; an electron has to receive energy equal to $E_g/2$ rather than E_g to cross the gap. Thus the study of the temperature dependence of electrical resistance allows us to find E_g . The most convenient way of doing this is to plot $\log(\rho)$ vs. $1/T$.

E.4.2 Metal–insulator transition

One type of a phase transition is that between metal and insulator or metal and semiconductor; it is most commonly called *Metal–Insulator Transition* [9, 10, 11] or MIT, even when the poorly-conducting phase is a semiconductor.

If we take a piece of metal, such as copper, and heat it up to its melting point, the melt will still be a metal. If we boil it (in vacuum) and evaporate the molten copper, the copper atoms constituting the vapor will be of course, strictly speaking, forming some sort of an insulator in a sense that they will not be conducting a current if placed in the electric field. This is not however what we call a metal-insulator transition. We are interested in a solid material which retains its general structure and composition, but, under the influence of an external parameter, for example, temperature undergoes an abrupt transition from a poorly-conducting semiconductor or insulator state to a well-conducting, metallic state or vice versa. There are some 50 or more compounds (but no pure elements) which exhibit MIT. For example, MIT was found in Fe_3O_4 , NiS , NbO_2 , V_2O_3 , etc. One such compound you will have a chance to study in this lab: it is vanadium dioxide, or VO_2 . All MIT's are first-order transitions, although this fact can be obscured by the presence of a disorder [9].

E.4.3 Essential physics of one type of MIT: the Mott transition

The theory of MITs is rather complex. What is more, several different MIT mechanisms may be at work in different compounds. Here we can only give a “taste” of a relevant theory by qualitatively describing the transition as envisioned by N. F. Mott [10, 11]. Other MIT models exist [9].



Figure 39: Sir Nevill Francis Mott (1905–1996), Nobel Prize of 1977.

Sir Nevill Mott, a famous British physicist, in 1949 discussed [10, 11] the following model: consider N hydrogen atoms are in a periodic lattice. If the distances between these atoms are small, N electrons which they contribute will fill half of the band which consists of $2N$ states (the reason there will be $2N$ states in such a band is that there are two states for each electron due to electron spin one-half: one up and one down [13, 15]). Thus such a hydrogen crystal will be a metal. Now Mott imagined increasing the distance between atoms (increasing lattice constant a). With increasing a , the band's energy span will decrease; it will become more and more narrow (consult solid state physics texts [13, 15] to see why), but, within the confines of the band theory, such “stretched out” hydrogen crystal will still be a metal with half-filled band for any, arbitrarily large, separation between atoms. This is of course nonsensical, because we know that eventually, at large enough atomic separations, the crystal must become an insulator.

This leads one to believe that, with increasing a , something new has to happen at a certain distance between atoms. Mott assumed that, starting from some critical value of a , electronic states no longer can be described as extended Bloch states. Electronic wavefunctions will become

localized. Such a crystal will be an insulator at $T = 0$: each electron will be localized at its respective atom. Such an insulator is called *Mott insulator*.

At sufficiently high T , however, such a material can conduct current, turning from an insulator into a metal. Indeed, consider two neighboring atomic sites, 1 and 2; each has an electron on it. If these two electrons have equally-oriented spins, the conduction is impossible in view of Pauli principle: an electron 1 cannot jump to site 2 which is occupied by electron 2 with the same quantum numbers. But suppose these electrons have opposite spins (what is called antiferromagnetic order). An electron from an atomic site 1 can in principle jump to the neighboring site 2, which has an empty energy state for electron 1's direction of spin (the other state being occupied by the electron with opposite spin which is sitting on site 2). If electron 1 jumps, it will leave behind a positively-charged hole on site 1. If such jumps are frequent, with an application of an electric field the current will flow, this current being carried by both electrons and above-mentioned holes they leave behind (holes moving in the direction of the field, and electrons against this direction, as usual [15]). To make this jump, however, an electron has to overcome Coulomb repulsion which exists between the two electrons residing on the same site — it has to have energy at least of the order of this Coulomb repulsion energy. In a Mott insulator, this Coulomb repulsion barrier energy serves as an analog of a forbidden gap E_g in a band picture of a dielectric. With raising temperature more electrons will be excited across this Coulomb gap, and thus the Fermi level (kinetic energy) will grow. Once the kinetic energy of electrons will become comparable to Coulomb repulsion energy, the current can flow. Thus, at a certain sufficiently-high temperature, Mott insulator will undergo a transition to a metallic state. Similar transitions can be envisioned as a function of doping, which likewise can increase electron concentration and thus kinetic Fermi energy, or pressure, which can reduce lattice constant a and, as a result, broaden the band and raise Fermi energy again.

Therefore, MIT in this picture takes place essentially as a function of free-electron concentration n , at a certain critical value n_C . External parameters which can drive n to this value of n_C can be temperature, pressure, doping, etc. Lattice transformations can also play a role too, leading to changes in n .

The nature of a MIT in our material, VO_2 , has been debated for many years. The latest opinion shared by the majority of workers in the field is that it is essentially a Mott transition as described above, and not a result of a lattice transformation. Instead, it is believed that lattice transformation itself is the secondary effect from an electronic Mott transition.

E.4.4 MIT in VO_2

The MIT in VO_2 was discovered in 1959, and extensively studied ever since. As a function of temperature it takes place at $T_C = 340 \text{ K} = 67^\circ \text{C}$ (transition is sharp in single crystals [17] but broad in polycrystalline samples [18]; you will measure it in thin films of VO_2 which exhibit a fairly broad transition). It is a first-order transition characterized by latent heat, density change, and lattice symmetry change, as well as by a dramatic change in resistivity. We have no way of detecting aspects of a transition other than resistivity in this lab.

You will measure resistance as a function of temperature. In Figure 40 we show an example of such a measurement taken from Ref. [18], performed on a 95 nm film of VO_2 . Your measurement will be different, as different films exhibit a somewhat different behavior; for example, the ratio of resistance just below the transition to resistance above the transition may be quite different in different samples. The overall features, however, are similar.

The transition, as is typical of all first-order transitions (see Section E.3.2), involves formation of a mixed-state, in which there are domains of metallic phase and domains of semiconducting phase. This is discussed in greater detail in Section E.4.5. The transition develops differently depending on the temperature increasing or decreasing — there is substantial hysteresis.

As the temperature increases from room temperature (25°C), resistivity will first decrease because our material is a semiconductor. But at some temperature, which is going to be close to

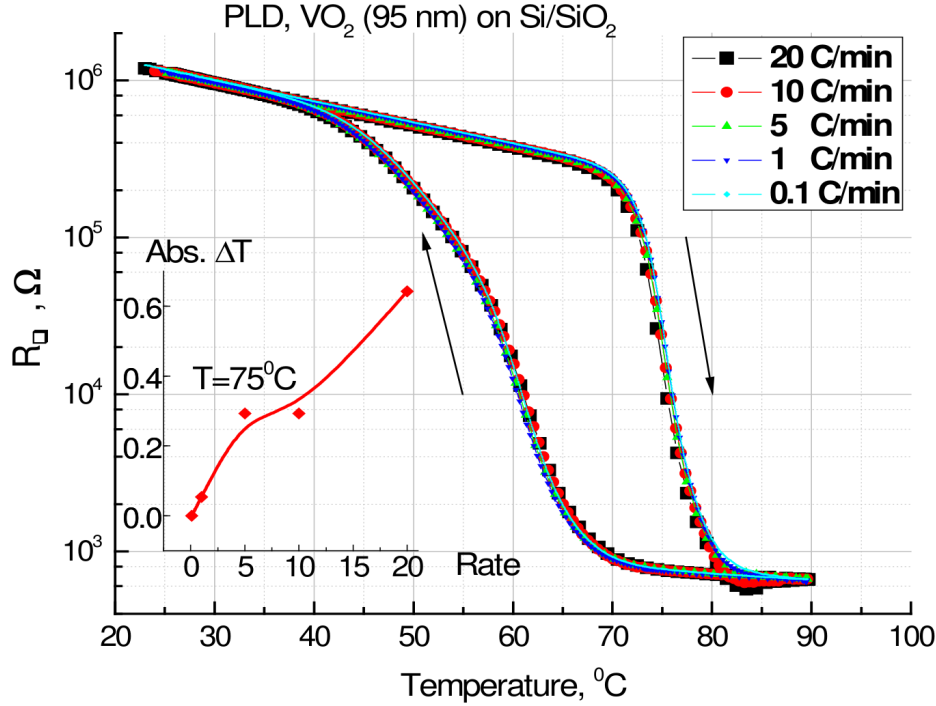


Figure 40: An example of a hysteresis loop in $R(T)$ measured on a 95 nm film of VO_2 deposited by pulsed laser deposition (PLD), from Ref. [18]. The figure shows several loops taken at different rates of temperature increase and decrease (temperature ramp rates); the upper inset gives a symbol key for the different ramps. The lower inset shows temperature shifts ΔT for different ramps. These are systematic errors in measuring temperature. Although on the logarithmic scale of the picture all these loops essentially overlap, there is a temperature shift of up to 0.6°C for the fastest ramp compared with the slowest. **Note that you will have to worry about such ramp rates and temperature shifts in your experiment.**

the $T_C = 67^\circ\text{C}$, metallic-phase domains will appear in an otherwise semiconducting material, and resistivity will start to decrease faster. Examination of Figure 40, which is plotted on a logarithmic scale in resistance, suggests that the first deviation from semiconducting behavior indeed takes place at about 67°C . The fraction of metallic domains grows with increasing T , and at some point a global metallic cluster connects the voltage leads on a sample; as a result all remaining semiconducting domains will be shorted out, and resistivity will drop to its metallic value. This is the result of the percolation transition [16], which, of course, is made possible by the appearance and growth of those metallic domains.

Question. If — for simplicity — we would imagine that resistivity of a semiconducting phase is infinitely greater than that of a metallic phase, where, at what temperature on the curve in Figure 40, would be the point at which a global metallic cluster is forming? This would be the point of a percolation transition. Read at least some of Ref. [16] and briefly discuss percolation physics in your report. At what fraction of metal phase is percolation transition expected to take place in a 2D (thin film) system?

If the temperature will be now lowered from about 90°C back to 25°C , the process will reverse, but with substantial hysteresis. All of this is shown in Figure 40. The transition, as can be seen, spans about 2.5 orders of magnitude in resistance. This measurement was performed on a VO_2 film; in single crystals resistance can change by 5 (!) orders of magnitude, and transition is quite sharp, with hysteresis of only about 1–2 degrees [17].

E.4.5 Development of a mixed state and hysteresis in VO₂

The text below is taken from Ref. [18], which has been slightly edited here. In this text you will find a qualitative explanation for the appearance of the coercive temperature T^* and of hysteresis, which was advanced by the authors of Ref. [18], primarily by Professor Serge Luryi (Chair, ECE Department, Stony Brook). Other explanations of the same phenomena were given by different authors, and they can be found in the literature [19].

The major issue in dispute here is whether each domain exhibits hysteresis on its own, or whether hysteresis is a collective property of the whole system, arising, as argued here, as a result of the boundary energy “cost”. Phases S and M refer to semiconducting and metallic phases. At the end there is a discussion of a phenomenon called *non-hysteretic branches*, which you can observe if you do appropriate measurements on your sample in this lab. Talk to the professor in charge of this experiment if you want to explore this phenomenon.

At a given temperature T inside the major hysteretic loop, some parts of the film have $T_C(x, y) < T$ and others $T_C(x, y) > T$. Here x, y are coordinates of a given location in a film. In the first approximation, the boundary wall between the S and M phases, at any given temperature T , is determined by the condition $T_C(x, y) = T$. In this approximation, the wall is highly irregular and its ruggedness corresponds to the scale at which one can define the local $T_C(x, y)$, i.e. to the characteristic length scale of the nanoscopic phase domains. On closer inspection, however, we need a refinement that takes into account the *boundary energy*, associated with the phase domain wall itself. The boundary energy is positive and to minimize its contribution to the free energy the domain walls are relatively smooth.

Let us examine the process of boundary motion. For concreteness, let us consider the heating branch. Below the percolation transition, M phase resembles lakes in the S phase mainland. With raising temperature the area of the M phase increases, lakes grow in size. When a boundary of a given lake is far from the other lakes, infinitesimal $\pm\delta T$ changes boundary length by infinitesimal amount, and the lake area by $\pm\delta A_M$. In other words, when the lakes are sufficiently separated, we envision a continuous, reversible, hysteresis-free process of $M \leftrightarrow S$ area redistribution, with neighboring configurations differing microscopically.

Let us now look at the formation of a link between two neighboring regions, which is the elementary step in the topological evolution of a global percolation picture. Let us focus on two metallic lakes that are about to merge. Since the boundary is smooth, at some temperature the distance between the lakes becomes smaller than the radius of curvature of either lake at the point they will eventually touch. Therefore, at some T the following two configurations will have equal energies: one comprising two disconnected M phase lakes that are near touching, but not quite, and the other with a finite link formed (Figures 41(a) and 41(b) respectively).

Both configurations are characterized by equal boundary lengths and therefore have equal free energy. The actual transition forming a local link, however, does not occur at that temperature because of an immense *kinetic* barrier (meaning: it takes significant atomic rearrangement) between these two *macroscopically* different configurations. The transition occurs at a higher $T_0 = T + T^*$, when it is actually forced, i.e. when the two phases touch at a point. Here T^* is the coercive temperature. In this picture the coercive temperature arises as a result of having a boundary between different phases; it does not pre-exist intrinsically within each domain. We associate the steep slopes of the major loop with the quasi-continuous formation of such links, i.e. with local topological changes, specifically with the merger of metallic lakes on the heating branch (HB) and of semiconductor lakes on the cooling branch (CB). Near the major loop ends T_S and T_M the global map consists of widely-separated M and S lakes, respectively. In these regions we expect to see non-hysteretic behavior, as was indeed observed.

Non-hysteretic branches (NHB). Consider now a small *backward* excursion from T_0 on the HB. As the temperature decreases, some of the M-phase recedes and the S-phase grows, changing the geometry of the global two-phase map. However, topologically, the last formed M-link does

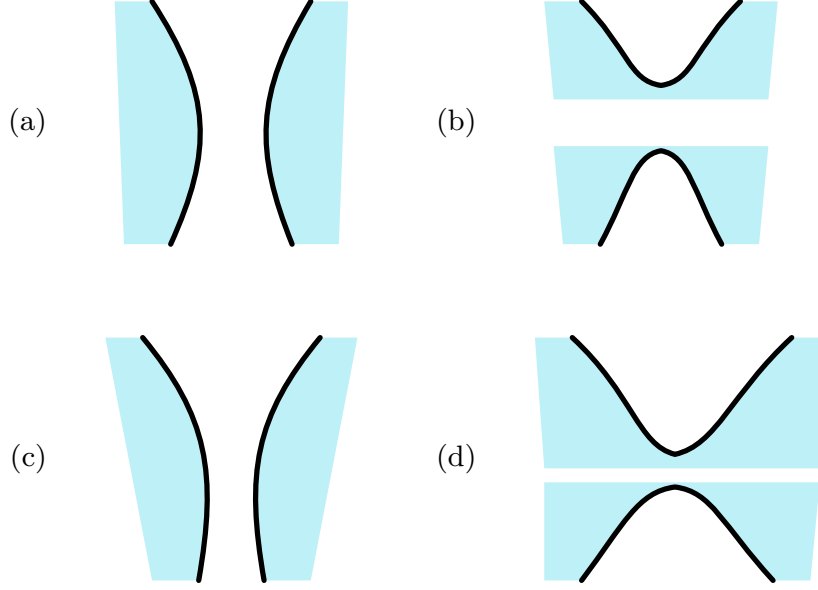


Figure 41: Semiconductor-metal boundary; the metallic phase is shaded. The top row, (a) and (b), corresponds to temperature T , and the bottom row, (c) and (d), to a higher temperature $T_0 = T + T^*$.

not disappear immediately for the same kinetic reason. One has two S regions that need to touch in order to wipe out the M-link. It takes a backward excursion of amplitude ΔT^* to establish an S-link and thus disconnect the last M-link. So long as we are within ΔT^* , i.e. stay on the same NHB, the area of S and M domains changes continuously, but the topology is stable and no new links are formed. Within the range of that stable or frozen topology, $\Delta T^* = \Delta T_{\text{NHB}}$, the resistivity of NHB will be single-valued and its T -dependence will be controlled by the percolating semiconductor phase. This explains the appearance of non-hysteretic branches and why they have semiconducting slopes.

E.5 BaTiO₃ and ferroelectricity

Ferroelectric materials [5, 6, 12, 13, 14] exhibit spontaneous electric polarization, the separation of the center of positive and negative electric charges. Recall the Landau theory conclusion discussed at the end of Section E.3.2, that the free energy will have a double minimum below T_C : in a ferroelectric state the electric polarization, which is the order parameter, is non-zero. This creates an excess positive charge on one side of the material, and negative charge on the other side. The polarization can be reversed by the application of an appropriate electric field: order parameter can be flipped from one minimum to the other with E-field. As we said above, ferroelectricity is named in analogy with ferromagnetism, which occurs in such materials as iron. The prefix *ferro*, meaning iron, was used to describe the property despite the fact that most ferroelectric materials do not contain iron.

In most ferroelectric materials the ferroelectricity disappears at high enough temperature, and the material turns into a regular dielectric. This state is called the *paraelectric*, in analogy to paramagnets. Most of the known transitions from the ferroelectric to the paraelectric state are second-order; the transition temperature is called the Curie temperature, in full analogy with the ferromagnetic transition. (Note: some ferroelectric materials melt before reaching this temperature.) In the paraelectric state these materials have very large dielectric constant: the temperature dependent dielectric constant $\varepsilon(T)$ increases as the temperature approaches the Curie temperature T_C , and shows a divergent behavior at $T = T_C$. We also note that near T_C

the dielectric constant is essentially equal to the susceptibility, a quantity you will find in many references on ferroelectricity. The relation between these two quantities is given in the next section.

Again, recall our — or rather, Landau’s theory’s — prediction made in Section E.3.2: we predicted very large polarizability near T_C , on both sides of a transition. Large polarizability means that induced charges in a dielectric will appear easily, compensating free charges on a capacitor plate — in other words, that relative dielectric constant k will be high.

Ferroelectricity is often accompanied by significant *piezoelectricity* (developing electric polarization under mechanical stress) and *pyroelectricity* (developing electric polarization when heated). Referring to our favorite (by now) Landau theory, we see that this is perfectly understandable: the system wants to fall into an energy minimum, and just about any outside stimulus will induce it to do so.

Ferroelectric materials have many applications, such as multilayer capacitors, gate dielectrics and waveguide modulators. Pyroelectricity is used in infrared detectors and other sensors. Piezoelectricity is widely applied in microphones, speakers and cigarette lighters, and in scientific applications such as the scanning tunneling microscope.

In this lab you will be studying a well-known ferroelectric material, an oxide of Ba and Ti, BaTiO_3 [12, 13, 20]. We note that oxides of this type are important and well studied in solid state physics and materials science, and not only because of their ferroelectric properties. These oxides all have a chemical formula of the type ABO_3 , where A is a metal with valency 1 or 2, either A^+ or A^{2+} ; B is a metal with valency 4 or 5, either B^{4+} or B^{5+} ; and O^{2-} is oxygen. The crystal structure of these materials is called the *perovskite* structure, after the naturally occurring mineral CaTiO_3 , which itself is called perovskite. Other examples of perovskite ferroelectrics include KNbO_3 , PbTiO_3 , KTaO_3 and $\text{PbZr}_x\text{Ti}_{1-x}\text{O}_3$ (the latter also called PZT).

When talking about crystal structure, we often show one chemical-formula unit of this structure, called the *unit cell*. The whole macroscopic (*bulk*, as it is often called) crystal of a material can be built by stacking unit cells in all three spatial dimensions.

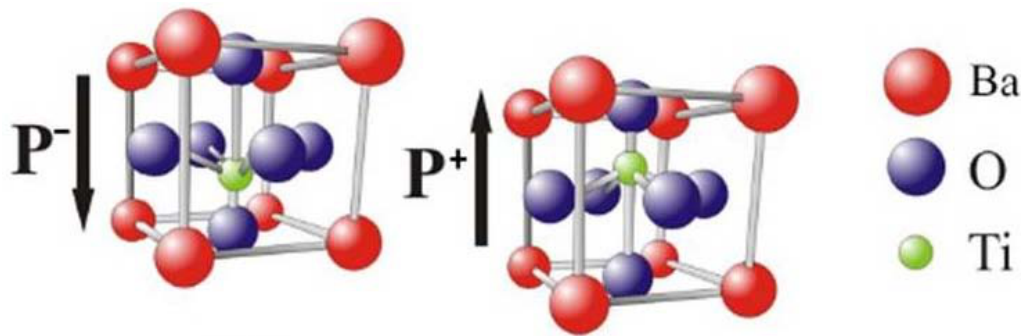


Figure 42: (a) Unit cell of a perovskite structure. (b) Ferroelectric distortion in BaTiO_3 .

Below a certain temperature T_C , some of the perovskites are characterized by the appearance of a *lattice distortion* which changes their structure from cubic, as in Figure 42(a), to a less symmetrical one (the structure which appears below T_C in our material is called *tetragonal*), as in Figure 42(b). Although all of the atoms constituting the unit cell move at least somewhat in this distortion, a rough picture of what happens can be obtained if you will visualize the central Ti^{4+} ion shifting within the unit cell by a small distance shown to be downwards in Figure 42(b); simultaneously, the four O^{2-} ions surrounding Ti^{4+} also shift, but in the opposite direction, i.e. upwards in Figure 42(b). All these shifts are quite small in comparison to the unit cell side a ; they are somewhat exaggerated in the picture above.

Note that the symmetry of the unit cell is reduced below the transition; a certain part of the total symmetry which existed above the transition is *broken* below the transition, as was discussed in Section E.3.2.

When the Ti ion is in the center of the cube, the charges are balanced, and the polarization is zero. If the Ti and oxygen ions move as illustrated in Figure 42(b) and in Figure 43, a finite *electric dipole moment* p appears in each of the unit cells, and the material becomes polarized; this polarization can be reversed with an application of an external field, as shown in Figure 43.

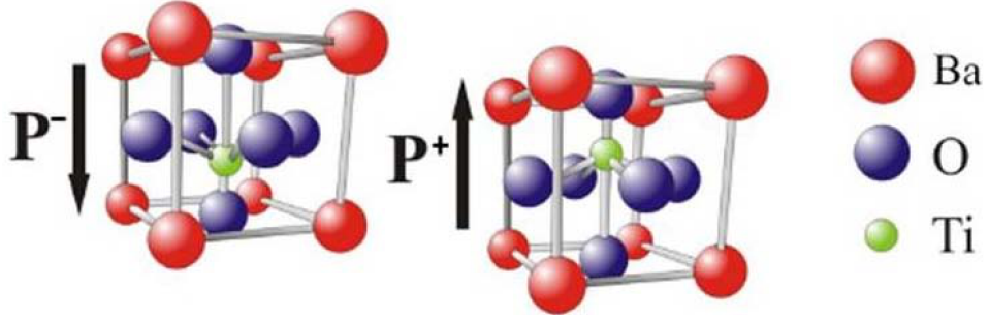


Figure 43: Unit cells of BaTiO₃ with the (green) Ti ion in two off-center positions: with p pointing up (called p^+ here) and with p pointing down (called p^- here). Below T_C , the dipole moments p^+ and p^- appear spontaneously without an external field. If there is an electric field E , for example pointing up, then p^+ will have lower energy than p^- .

In this lab we will use a temperature-controlled environment to study the dielectric constant as the sample undergoes the phase transition, and we will search for signs of “divergence” in the dielectric constant. We will also attempt to extract some microscopic information from these measurements, as described in the main text of this write-up and also below.

A note on terminology and notation. The dielectric properties of any material depend on the frequency at which the properties are measured. This is described by the dielectric function $\varepsilon(\omega)$. In our case the properties depend on the temperature as well, and therefore we have a two-variable function $\varepsilon(\omega, T)$. We use the term “temperature-dependent dielectric constant” to represent the low-frequency limit of the dielectric function, $\varepsilon(T) = \varepsilon(\omega = 0, T)$. The frequencies applied in our studies (25 kHz) are actually in this low frequency limit. Furthermore, we will distinguish between dielectric function $\varepsilon(T)$ and dielectric constant $k(T)$, which is just the relative dielectric function $k = \varepsilon(T)/\varepsilon_0$. We may call it a “constant”, following the normal usage, but in our case, when studying a ferroelectric, we will recognize that it is a temperature-dependent quantity.

E.6 Paraelectric to ferroelectric transition; Curie–Weiss law

BaTiO₃ undergoes a second-order paraelectric to ferroelectric phase transition upon cooling through $T_C = 120^\circ\text{C}$. Above that temperature it has cubic structure and no local dipole moment; below it transforms to a new, slightly-non-cubic, tetragonal structure in which Ti ion and oxygen ions move in opposite directions creating a local dipole moment in each unit cell. In other words, each unit cell becomes electrically polarized and, below T_C , these local dipoles align to form a macroscopic polarization P .

As was stated in Section E.3.2, electric susceptibility is defined as $\chi_E = dP/dE$, or, when P is linear with E , as $\chi_E = P/E$. We mentioned in Section E.3.2 that Landau theory predicts that in the vicinity of a phase transition χ_E is large (in fact, diverging at $T = T_C$) and that it follows Curie–Weiss Law

$$\chi_E = \frac{C_{CW}}{T - T_C}. \quad (28)$$

The relationship between χ_E and the dielectric constant ε is given by $\varepsilon = \chi_E + 1$ [2, 3, 6, 13]. Near the transition χ_E is large, and so, to a very good approximation, $\varepsilon = \chi_E$. Therefore

Curie–Weiss law (28) can be re-written for the relative dielectric constant as

$$k = \frac{C_{CW}}{\varepsilon_0 (T - T_C)}. \quad (29)$$

Again, to clarify the meaning of symbols and to avoid confusion, in this equation k is a *relative* dielectric constant, in other words, it is equal to the ratio of capacitance with a dielectric to the capacitance with the same geometry but with vacuum between the plates, $k = C/C_0$. As was mentioned above, although it is called a “constant”, it is not a constant with respect to temperature. The formula for a parallel-plate capacitance with vacuum (or air) between the plates is $C_0 = \varepsilon_0 A/d$, where A is one plate’s area, and d is the distance between the plates. The formula for the capacitance with dielectric is $C = \varepsilon A/d = \varepsilon_0 k A/d$. Here ε_0 , the permittivity of free space, is equal to $8.85 \times 10^{-12} \text{ C}^2 \text{ N}^{-1} \text{ m}^{-1} = 8.85 \times 10^{-12} \text{ F m}^{-1}$.

In Section E.3.2 we saw that Landau theory predicts this law both above and below the transition: we had, from (25), $\chi_E = C_{CW1}/(T_C - T)$ for $T < T_C$ and $\chi_E = C_{CW2}/(T - T_C)$ for $T > T_C$, with $C_{CW1} = 1/(2a)$ and $C_{CW2} = 1/(4a)$, a being expansion coefficient. This law is indeed observed in the vicinity of ferromagnetic and ferroelectric transitions, where experimental data mostly concentrates on paramagnetic or paraelectric states at $T > T_C$. You will have a chance to see it in your own measurements, and even to test it on both sides of the transition.

Problems. As was already discussed in the main text of this write-up, your sample may have a somewhat irregular shape, so you will at best be able to estimate k perhaps within a factor of 2. Additionally, you are measuring a ceramic sample, which will tend to underestimate intrinsic k in your material, as was discussed in Section 9. The latter problem can be avoided or corrected, at least in part, by making a correction as described in that section (let us call this a *capacitance correction*). You want to use your measured C values to calculate the value of C_{BTO} , which is the intrinsic capacitance of your material. Thus when we will be talking about your capacitance data, we will assume that you made the suggested capacitance correction. The k values which you will be using in the k vs. T data should be using this correction.

Keeping in mind these caveats, by fitting your corrected k vs. T data to formula (29) you will be able to find (or at least to estimate) the value of C_{CW} . You may, for example, plot $1/k$ vs. $T - T_C$ and find the slope of that line, which is $1/C_{CW}$. You will be using your measured T_C here, which is easy to identify as the temperature of the sharp cusp in your data.

Ramp slowly. You want to make sure that your temperature ramp is slow enough to trust your T_C measurement.

The temperature interval of your data above T_C will be rather narrow, but you will be able to estimate C_{CW} . It will be interesting to also fit the data below the transition, to see if you can verify the theoretical factor of two difference between C_{CW1} and C_{CW2} (see Section E.3.2). As we said above, the largest uncertainties in the value of C_{CW} will come from not knowing the capacitor geometry well, and from measuring a ceramic sample (although the latter may be in part remedied by the capacitance correction), rather than from your fits to a Curie–Weiss law.

Why do we want to find C_{CW} ? Because there is a theoretical connection between the value of C_{CW} and microscopic parameters, and we will try to obtain it in the next section.

E.7 Derivation of C_{CW} in terms of a microscopic dipole moment

In a ferromagnetic material microscopic magnetic moments μ (fundamentally, spins) are always present; the difference between paramagnetic and ferromagnetic state is merely in the organization, alignment or misalignment of these pre-existing moments. In a paraelectric material, a unit cell is not carrying its p in the same way as an atom carries its spin; in fact the microscopic electric dipole moment p above the T_C should not exist. The lattice symmetry changes at T_C , and above T_C we have a cubic cell with Ti in the center, so that microscopic $p = 0$.

Here we wish to make the following comment: Landau's potential above T_C (Figure 37(a)) has only one minimum, and this may suggest that all electrical dipoles p and magnetic moments μ must disappear above T_C . Yet we know that individual spins do not disappear at $T > T_C$ in a paramagnetic state. Thus we realize that the notion of order parameter m (in our case, polarization P) is macroscopic in nature; it does not tell us anything about the existence or non-existence of microscopic p 's above the T_C . Our conclusion about local $p = 0$ in a paraelectric state in BaTiO_3 is based on the known structure and symmetry of its unit cell [12, 13] and not on Landau's theory.

Another difference between magnetic and electric cases is in that quantum mechanics imposes certain requirements on the quantum numbers and energy levels in the magnetic case: as is well-known, a degenerate single energy level in $B = 0$ field will split into two levels for spin one-half electrons in $B \neq 0$ field [6, 7]. This may not be the same in case of p 's and E -field: there is no quantum number associated with p . Yet we will also use the same two-level picture here. In other words, we will assume that microscopic dipole p can point either along the E field or against it, as shown in Figure 43, and will ignore its other possible orientations with respect to the direction of E .

The analysis of the ferromagnetic case can be found in a number of books, e.g. in Refs. [5, 6, 13, 15]. Our derivation here is similar to these treatments, but somewhat simplified.

Let us therefore consider a material in a ferroelectric state, at $T < T_C$, where we can assume that one dipole moment p exists in each unit cell of volume a^3 ; thus in 1 m^3 there are $N = 1/a^3$ microscopic dipoles p . Let us roughly estimate the value of a microscopic dipole moment in BaTiO_3 . The "generic" dipole is shown in Figure 44.

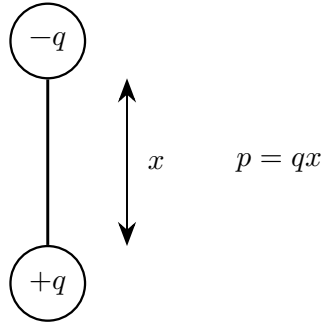


Figure 44: A generic electric dipole, $p = qx$.

It consists of two opposite electric charges, $-q$ and $+q$, separated by a distance x ; the formula is $p = qx$. Whatever ionic displacements happen inside the unit cell of a ferroelectric, and whatever charge imbalance they will create locally, the total charge cannot change, and the unit cell must remain neutral. Thus our effective electrical dipole which represents this polarization, must be neutral as well, with $-q$ and $+q$ adding to zero. In this formula we treat q as some effective electric charge, which, we recognize, must be of the order of electron charge e , or, at most, it can be several times larger than e . Indeed, we know that in BaTiO_3 , Ti ion carries charge $+4$, and Ba ion carries charge $+2$, while three oxygen ions carry charge -6 , with the sum of all these charges being, of course, zero. We will therefore take $q = 6e$ in our effective dipole.

The second quantity in the dipole formula is x , an effective displacement which represents charge-weighted displacements of all of the charges in the unit cell. In BaTiO_3 , in the first, rough approximation, we can think of x as a displacement of a Ti^{4+} ion, which in a cubic state resides in the middle of a cell and shifts up or down (Figure 43) in a ferroelectric state. Therefore, the absolute upper limit on the value of x is $a/2$, half size of the unit cell. In reality we expect x to be much smaller than this upper limit. The lower limit is hard to predict; it can be a small fraction of a . So we expect to find $0 < x < a/2$. The length of a cubic unit cell side in BaTiO_3 is

$a = 0.4 \text{ nm}$ [12, 13, 14], which gives us $0 < x < 0.2 \text{ nm}$.

Thus we have a rough upper-limit estimate

$$p_{\max} = q x_{\max} = 6e x_{\max} \approx (9.6 \times 10^{-19}) \times (0.2 \times 10^{-9}) \approx 2 \times 10^{-28} \text{ C m}. \quad (30)$$

We also know that in the electric field E , potential energy of p is $U = -(\mathbf{p} \cdot \mathbf{E})$, i.e. it changes from the minimum value ($-pE$) for p being parallel to E to the maximum value ($+pE$) for p being anti-parallel to E . The thermal energy at room temperature is $k_B T = 1.38 \times 10^{-23} \times 300 = 4.1 \times 10^{-21} \text{ J}$. We see that it will take electric field of $2 \times 10^7 \text{ V m}^{-1}$ to make maximum dipole energy equal to thermal energy. We conclude that for any reasonable value of electric field E , $U \ll k_B T$.

Now, near the transition, where macroscopic polarization P (the order parameter) is small, in the absence of an external field E , all microscopic p 's have the same (zero) energy. Once the external field E is applied, this energy will split into two levels, one above and one below the initial level, the energy difference between these two levels being $\Delta U = pE - (-pE) = 2pE$, as shown in Figure 45.

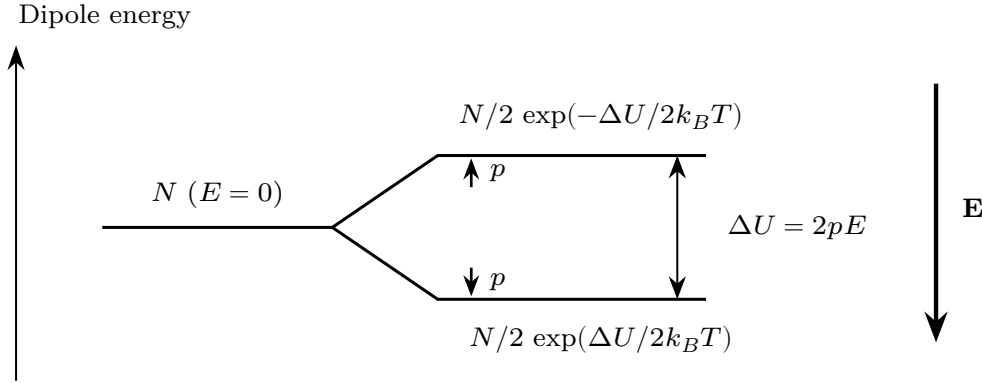


Figure 45: Splitting of the dipole energy level in an applied electric field.

We still have $\Delta U \ll k_B T$, which means that there will be approximately $N/2$ dipoles on each of the two energy levels. Nonetheless, the small ΔU is responsible for the population difference between the two levels: the lower level will have

$$N_1 = \frac{N}{2} \exp\left(\frac{\Delta U}{2k_B T}\right) \approx \frac{N}{2} \left(1 + \frac{\Delta U}{2k_B T}\right) \quad (31)$$

dipoles, while the upper level will have

$$N_2 = \frac{N}{2} \exp\left(-\frac{\Delta U}{2k_B T}\right) \approx \frac{N}{2} \left(1 - \frac{\Delta U}{2k_B T}\right) \quad (32)$$

dipoles, so that $N_1 > N_2$. The total number of dipoles is of course still $N_1 + N_2 = N$. The relatively small difference $N_1 - N_2$ is what is responsible for macroscopic polarization P . So we have

$$N_1 - N_2 = \frac{N}{2} \frac{\Delta U}{k_B T} = \frac{N}{2} \frac{2pE}{k_B T} = \frac{NpE}{k_B T}, \quad (33)$$

Macroscopic polarization is

$$P = p(N_1 - N_2) = \frac{Np^2 E}{k_B T}. \quad (34)$$

Finally we will get for susceptibility $\chi_E = P/E$, and for relative dielectric constant,

$$\chi_E = \frac{Np^2/k_B}{T}, \quad k = \frac{Np^2/\epsilon_0 k_B}{T}. \quad (35)$$

The one difference between the last formula and formula (29) is in the temperature dependence: while in (29) we had $1/(T - T_C)$ (or, below the transition, $T_C - T$), here we got simply $1/T$. But we did not develop in this section a full theory of a phase transition; we just assumed that transition takes place at some temperature; in fact, we could assume that it took place at zero temperature; then $T - T_C$ would be the same as T_C . So we should not be very surprised that we did not get $1/(T - T_C)$; the correct temperature dependence was addressed in Landau theory in Section E.3.2. You can also find a derivation of these formulas with the correct $1/(T - T_C)$ dependence in Ref. [6].

So we see that we should identify the numerator with C_{CW} ; by doing so we obtain

$$C_{CW} = \frac{Np^2}{\varepsilon_0 k_B}, \quad (36)$$

or, in terms of a microscopic dipole moment p ,

$$p = \sqrt{\frac{\varepsilon_0 k_B C_{CW}}{N}}, \quad (37)$$

where $N = 1/a^3$ and C_{CW} is a measured quantity which can be obtained from fitting your k vs. T data to the Curie–Weiss Law.

As an exercise, check the dimensionality of C_{CW} in (36). It should come out to be K (kelvin), so that k in formula (29) will indeed be a *relative* quantity without units.

What you should calculate. You should find C_{CW} derived from your data (again, we assume that you applied the capacitance correction), substitute it in formula (37), and calculate p , and then compare it to our upper-limit estimate $p_{\max}$ for this quantity. Once you will know p , you can calculate x , the effective dipole displacement. This last value can be compared to the literature: indeed, transition-induced displacements of various atoms in the unit cell have been studied by X-rays. You should also find in the literature the value of C_{CW} for single-crystalline BaTiO_3 , and, using it in formula (37), calculate again p and x . This time the values will be more realistic, as you will be using the correct, intrinsic value of C_{CW} . Compare the results with your own, comment on discrepancies.

References

- [1] *Model SR830 DSP Lock-In Amplifier*, Stanford Research Systems, revision 2.5 (2011). A copy is in the PyFERRO repository, pyferro/docs/manuals.
- [2] *Lock-in amplifier*, Wikipedia. http://en.wikipedia.org/wiki/Lock-in_amplifier — about lock-in amplifiers; at the end of this article you will find references to more lock-in tutorials.
- [3] *AC circuits*, School of Physics, University of New South Wales. <http://www.animations.physics.unsw.edu.au/jw/AC.html> — a tutorial on AC circuits, with animations.
- [4] *PID controller*, Wikipedia. http://en.wikipedia.org/wiki/PID_controller — about temperature control.
- [5] M. Gitterman and V. Halpern, *Phase Transitions: A Brief Account with Modern Applications*, World Scientific (2004).
- [6] P. Papon, J. Leblond and P. H. E. Meijer, *The Physics of Phase Transitions*, Springer (2001).
- [7] *Landau theory of phase transitions*, lecture notes, University of Helsinki. <http://theory.physics.helsinki.fi/~stafyi/SPI06chap10.pdf>
- [8] *Second order phase transitions*, P. Hadley, TU Graz. http://lamp.tu-graz.ac.at/~hadley/ss2/landau/second_order.php
- [9] V. Dobrosavljević, *Introduction to Metal–Insulator Transitions* (2011), arXiv:1112.6166v1.
- [10] N. F. Mott, *Metal–Insulator Transitions*, Taylor & Francis Ltd. (1974), ISBN 0-85066-079-3.
- [11] N. F. Mott, *Proc. Phys. Soc. (London)* **A62**, 416 (1949).
- [12] F. Jona and G. Shirane, *Ferroelectric Crystals*, Pergamon Press / The Macmillan Co. (1962). A copy of the relevant pages is provided in the lab.
- [13] J. S. Blakemore, *Solid State Physics*, 2nd edition, W. B. Saunders Co. (1974), Chapter 5, “Dielectric and magnetic properties of solids”.
- [14] *Phase transitions in ferroelectrics*, Advanced Lab, University of Michigan. http://instructor.physics.lsa.umich.edu/adv-labs/Ferroelectrics/phase_transitions.pdf
- [15] N. W. Ashcroft and N. D. Mermin, *Solid State Physics*, Thomson Learning, Toronto (1976).
- [16] *Percolation*, Chapter 13 of course notes, Clark University. <http://physics.clarku.edu/courses/125/gtcdraft/chap13.pdf>
- [17] C. N. Berglund and H. J. Guggenheim, “Electronic properties of VO₂ near the semiconductor–metal transition”, *Phys. Rev.* **185**(3), 1022–1033 (1969).
- [18] M. Gurvitch, S. Luryi, A. Polyakov and A. Shabalov, “Nonhysteretic phenomena in the metal–semiconductor phase-transition loop of VO₂ films for bolometric sensor applications”, *IEEE Transactions on Nanotechnology* **9**, 647–652 (2010).
- [19] For example, see T. G. Lanskaya, I. A. Merkulov and F. A. Chudnovskii, “Hysteresis effects at the semiconductor–metal phase transition in vanadium oxides”, *Sov. Phys. Solid State* **20**(2), 193 (1978), and a number of more recent papers treating hysteresis in VO₂.

- [20] *Ferroelectric materials: phase changes*, DoITPoMS, University of Cambridge. http://www.doitpoms.ac.uk/tlplib/ferroelectrics/phase_changes.php — about the structure of BaTiO_3 , with very nice graphics and illustrations.

Sections relating to ferroelectricity in BaTiO_3 were created by Laszlo Mihaly and in part by Michael Gurvitch; sections relating to phase-transition theory and the metal–insulator transition in VO_2 were created by Michael Gurvitch, both in 2012. Thanks to Xu Du, Matthew Dawber, Alexander Shabalov and Carlos Marques for contributions to the design and realization of the measurement and for useful comments. The 2026 edition of this manual was typeset into $\text{\LaTeX}$ with Claude Opus 4.5 and updated by Logan Levack and Qiyun Ma.